

BENOIT SYSTEMES

LIGHT ASSIST 3 LIGHT ASSIST 3 XL

USER MANUAL

September 2025






Follow us!

1 SETTING UP THE POWER ADD-ON

Your wheelchair has been *fitted with hook bars* and hooks* that will allow you to set up your power add-on very quickly.


1. Put the brakes on the wheelchair to make the setting up easier.
2. Lift the power add-on by the upper handle and place the axles on the hooks.
3. Release the handle and check that the power add-on is locked onto the hooks.

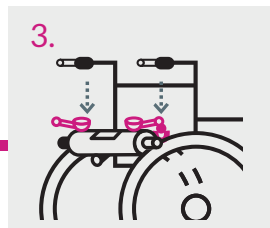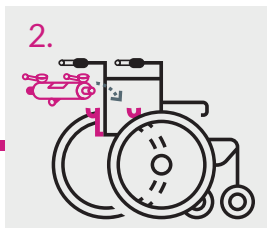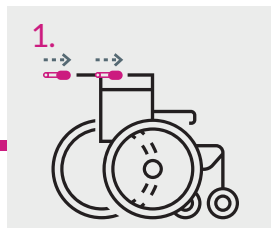

**removable for comfort chairs*

2 TURN THE HANDLES TO ENGAGE


1. Extend the telescopic handles to gain more force.
2. Engage both sides.

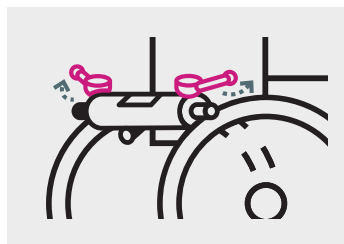

EXTRA FEATURE

Simple shift of the clutch handles is enough to disengage the power-add-on and get your manual wheelchair back.

SWITCH FROM POWER TO MANUAL MODE

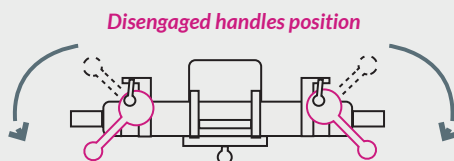


Telescopic handles


- 1 Extend the handles.
- 2 Turn the handles until they reach the stop.

3 SETTING UP THE PUSH HANDLES

Brackets have been fitted on the wheelchair rear handles to put the Light Assist handles. Two or three* handles are supplied with the power add-on:

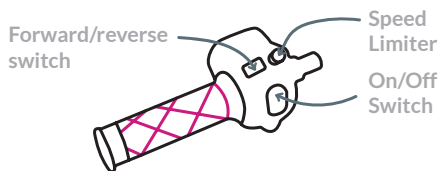

A throttle handle to drive the chair.

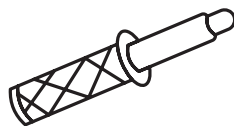



One or two static handles, to manually push the wheelchair when needed.

Handles are locked with a pin at the left or right-hand side on the wheelchair rear handles.

**for comfort chairs*

4 SETTING UP THE BATTERY


1. Place the battery on the battery sliders on the power add-on.
2. Firmly push the battery until you hear the connectors slot into each other.

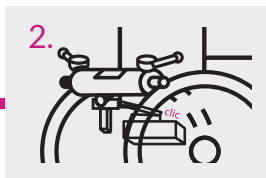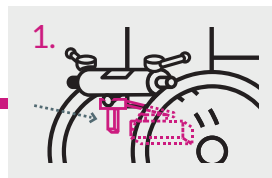

5 SETTING UP THE ANTI-TIP

1. Clip the anti-tip into the rear handle.
2. Check the mounting orientation (wheel bottom facing) and the pin locking.

Light Drive is now ready for use!

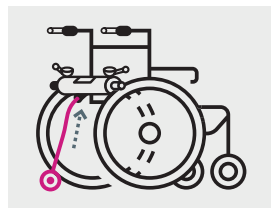

INSTRUCTION TIPS

Here are a few tips to ensure a correct use and a long life of your power add-on:



Charge your battery after each day of use.

- ⓘ The battery must be switched on to charge*.

***For more details, refer to the battery user manual.**



Inflate your tyres regularly. The operating pressure should be between 6 and 8 bar.






Regularly clean the rollers of the power add-on.

Follow our [Youtube channel](#) to get regular advice.

6 REMOVING THE POWER ADD-ON


1. Put the brakes on the wheelchair.
2. Keep the handles engaged.
3. Remove the anti-tip.
4. Hold the battery handle and pull it backwards.
5. Remove the control handles and place them in the storage holders on the power add-on.
6. Disengage the power add-on by turning the handles.
7. Lift the motor using the upper handle.

 The device can only be removed by lifting the upper handle.

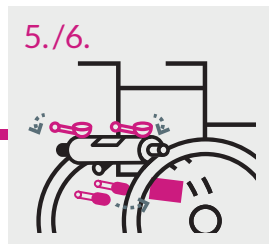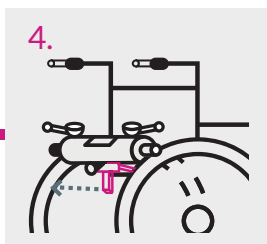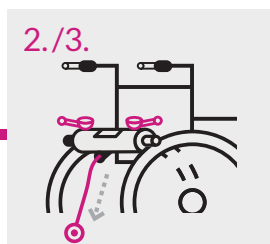


7 CHARGING THE BATTERY

If the battery remains on the wheelchair, connect the charger to the connector on the back of Light Assist.

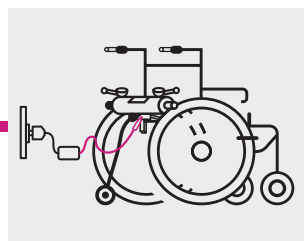

If the battery is removed, connect the charger to the red connector on the battery.

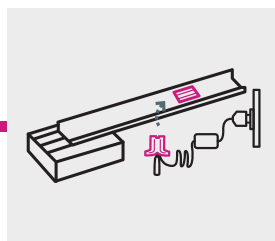

**READ THE
SERVICING
MANUAL AND
BATTERY USER
MANUAL**

for driving, security
and maintenance
instructions

**Then plug the charger into the power supply.
Disconnect once the load is complete.**

This medical device is a regulated health product which, under these regulations, bears the CE mark. In its capacity as manufacturer, BENOIT SYSTEMES declares that this product is a Class I device, compliant with the safety and performance requirements of European Regulation 2017/745/EEC on medical devices.

📍 7 Rue du Pont, 21450 Billy-Lès-Chanceaux FRANCE
☎ (+33) 3 80 96 51 25 | @ contact@benoitssystemes.com

WWW.BENOITSYSTEMES.COM