



LIGHT DRIVE 3

LIGHT DRIVE MINI 3 / LIGHT DRIVE 3 XL

USER MANUAL

January 2026






Follow us!

1 SETTING UP THE POWER ADD-ON


Your wheelchair has been fitted *with hook bars and hooks* that will allow you to set up your power add-on very quickly.

1. Put the brakes on the wheelchair to make the setting up easier.
2. Lift the power add-on by the upper handle and place the axles on the hooks.
3. Release the handle and check that the power add-on is locked onto the hooks.

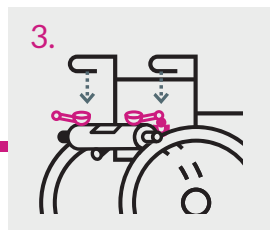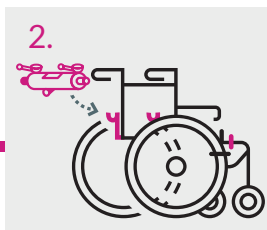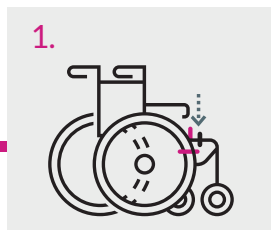


2 LOCKING THE HANDLES

1. Extend the telescopic handles to gain more force.
2. Engage both sides.

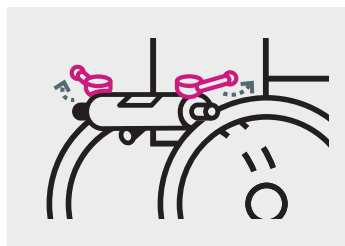

Simple shift of the clutch handles is enough to disengage the power add-on and get your manual wheelchair back.

SWITCH FROM POWER TO MANUAL MODE

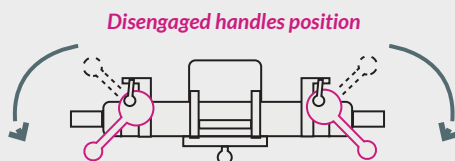


- ① *Extend the handles.*
- ② *Turn the handles until they reach the stop.*

3 PLACING THE JOYSTICK

A **controller bracket** at the front of the wheelchair or under the armrest will hold the joystick.




1. Place the joystick in the controller bracket.
2. Place the joystick cable alongside the seat to avoid friction on the wheel.

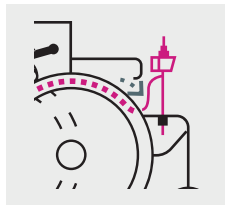

To retract the joystick, press or pull the button on the bracket. Lift it up until the stop to get back to original position.

4 SETTING UP THE BATTERY

1. Check that the battery is switched on*, then place it on the battery sliders on the power add-on.
2. Firmly push the battery until you hear the connectors slot into each other.

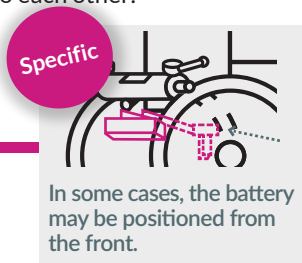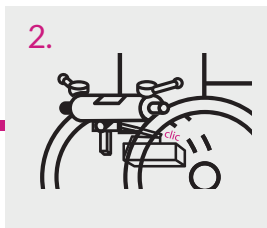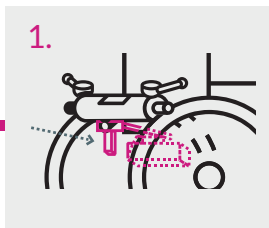


In some cases, the battery may be positioned from the front.

5 SETTING UP THE ANTI-TIP

1. Clip the anti-tip into the rear handle.
2. Check the mounting orientation (wheel bottom facing) and the pin locking.

Light Drive is now ready for use!

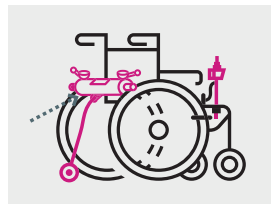

INSTRUCTION TIPS

Here are a few tips to ensure a correct use and a long life of your power add-on:



Charge your battery after each day of use.

- ⓘ The battery must be switched on to charge*.

*For more details, refer to the battery user manual.



Inflate your tyres regularly. The operating pressure should be between 6 and 8 bar.






Regularly clean the rollers of the power add-on.

Follow our **YouTube channel** to get regular advice.

6 REMOVING THE POWER ADD-ON


1. Put the brakes on the wheelchair.
2. Keep the handles engaged.
3. Remove the anti-tip.
4. Hold the battery handle and pull it backwards.
5. Remove the joystick and place it in the storage bracket on the power add-on.
6. Disengage the power add-on by turning the handles.
7. Lift the motor using the upper handle.

 The device can only be removed by lifting the upper handle.

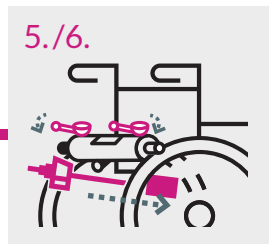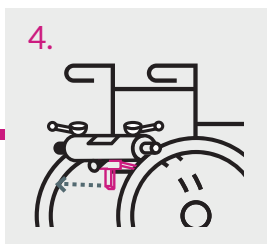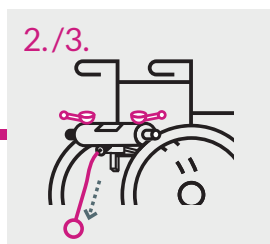


7 CHARGING THE BATTERY

If the battery remains on the wheelchair, connect the charger to the joystick.

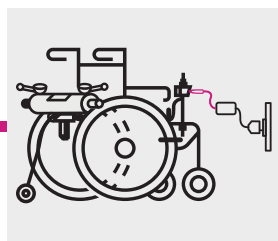

If the battery is removed, connect the charger to the red connector on the battery.

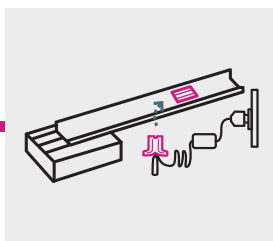

**READ THE
SERVICING
MANUAL AND
BATTERY USER
MANUAL**
for driving, security
and maintenance
instructions

**Then plug the charger into the power supply.
Disconnect once the load is complete.**

This medical device is a regulated health product which, under these regulations, bears the CE mark. In its capacity as manufacturer, BENOIT SYSTEMES declares that this product is a Class I device, compliant with the safety and performance requirements of European Regulation 2017/745/EEC on medical devices.

📍 7 Rue du Pont, 21450 Billy-lès-Chanceaux FRANCE
☎ (+33) 3 80 96 51 25 | @ export@benoitssystemes.com

WWW.BENOITSYSTEMES.COM