



DICHIARAZIONE DI CONFORMITÀ UE / UE CONFORMITY DECLARATION



PRODOTTO / PRODUCT	JOB
IDENTIFICATIVO ATTRIBUITO DAL FABBRICANTE / MANUFACTURER PART NUMBER	REF S105
EMDN	Y122103 PUSH WHEELCHAIRS
UTILIZZO PREVISTO / INTENDED USE	Il prodotto è destinato a essere impiegato per attenuazione o compensazione di una lesione o di una disabilità. In particolare la carrozzina è intesa per essere utilizzata da persone con difficoltà di movimenti temporanee o permanenti e che hanno la necessità di compiere spostamenti sulla spiaggia e su altri tipi di superfici impervie. The product is intended to be used for alleviation or compensation for an injury or disability. Particularly the wheelchair is intended to be used by those groups of users with temporary or permanent mobility difficulties confined to a sitting position who need to move on the beach and other types of impervious surfaces.
FABBRICANTE / MANUFACTURER	Neatech.it s.r.l. via Antonio de Curtis 4/A – 80040 Cercola (NA) – Italia Tel. +39 081 5551946 – info@neatech.it – www.neatech.it P. IVA IT04812481218 – REA NA 715393 Eudamed actor ID/SRN IT-MF-000031601
UDI DI BASE / BASIC UDI DI NUMBER	805538493S105CC
UDI DI / UDI DI NUMBER	08055384930025



We have great concerns about quality, safety and usability of our products. We are always working to improve them. Please download last version of this manual at this link.

<https://www.neatech-download.it/s105-manual/>

Index

1	PRODUCT INFORMATION	1
2	PREPARATION FOR FIRST USE	2
2.1	Mounting the wheelchair	2
3	ADJUSTMENTS	6
3.1	Legrest adjustment	6
4	USE OF THE PRODUCT	7
4.1	Armrests.....	8
4.2	Third wheel kit	8
4.3	Sun awning.....	10
4.4	Other accessories	13
4.5	Transport and storage	13
5	MAINTENANCE	14
5.1	Cleaning	14
5.2	Controls to be performed on the product.....	15
5.3	Tire puncture	16
5.4	Reuse	16
6	SPECIFICATIONS	17
6.1	Dimensions	17
7	WARRANTY TERMS	18
7.1	Serial number.....	19

Symbols in this manual



WARNING

This symbol means presence of danger for the user or damage for the product. Always follow instructions when this symbol is present.



INFORMATION

This symbol means general information intended to simplify or best explain the use of the product.



CONTACT INFORMATION

This symbol means the need of contacting an authorized service center or the manufacturer.

	INCIDENT REPORTING <u>If an incident occurs, please contact the manufacturer.</u>
	In order to obtain product safety notices and product recalls, for example by ensuring the supplier has up-to-date contact details visit the website of the manufacturer.
	AUTHORIZED SERVICE CENTER In order to obtain a list of authorized service center, please contact the manufacturer.
	This product and all its components can not be treated as household waste. For more detailed information on how recycling and disposal this product contact your local waste disposal service.
	Any unauthorized modification to the product may increase the risk of personal injury and damage to the product itself. All modifications should be done by an authorized service center. Do not use any unauthorized accessories or spare parts on the product.
	Do not use the product in combination with other medical devices without first having considered any risk due to combination of more products. Please pay special attention to the combination of flame-resistant and non-flame-resistant materials.
	<u>WARNING: It is prohibited to use the product or its parts for any purpose other than that indicated. For a correct use please follow the instructions given in this manual. The manufacturer disclaims any responsibility for damages caused by improper use of the product.</u>
	Information in this manual may be subject to change without notice. All information, pictures and specifications are based on the product details that were available at the time of preparation of this document. They are representative examples and they are not intended to be exactly as the actual product.
	TIPPING HAZARD Tipping hazard is strongly reduced because of the design of the product according to EN 12182. In any case, please pay special attention during the adjustments and use of the product to prevent any damage to the user or product itself.
	Center of balance of the device and so its stability can be affected by user position.



TEMPERATURE

The temperature of some surfaces may increase when the product is exposed to external heat sources as direct sunlight.



PINCH HAZARD

Make sure you have proper space in areas you will travel through to minimize pinching or entrapment of body parts.



INFORMATION

For information on how to obtain information and instructions in a format appropriate for use by visually impaired people please contact the manufacturer.



Do not operate the product if it is behaving abnormally or erratically.



Do not stand on the product. Always use caution when transferring in or out of the seat. Every precaution should be taken to reduce the transfer distance. Do not let children use the device without supervision.



Do not install, maintain or operate the product without reading all warnings and this entire manual.

Always keep this manual in connection with the product.



Do not carry passengers on the product independently of the age of the passenger. The product is not designed for weight training and is unsafe for use as a seat while weight training. Do not shift your weight or sitting position toward the direction you are reaching as the product may tip over. Do not stand on the frame of the product. **Some pathologies may limit your ability to use your device safely. Be sure to consult with a doctor about your physical limitations.**



PELVIC BELT

The wheelchair has the predisposition for a pelvic belt. **Pelvic belt is only design to position the user and not for any protection in case of accident.**



Dealing with uphill

The product is not explicitly designed to deal with uphill. Please pay special attention when dealing with uphill.



Dealing with downhill

The product is not explicitly designed to deal with downhill. Please pay special attention when dealing with downhill.



Dealing with side slope

The product is not explicitly designed to deal with side slope. Please pay special attention when dealing with side slope.



Turning with the product

When turning with the product, always use great caution.



Dealing with obstacles

The product is not explicitly designed to deal with obstacles. Please pay special attention when dealing with obstacles.



WARNING

The stopping distance on slopes can be significantly greater than on level ground



Lift the product

Do not lift the product with the user on board. If you really need to lift the product, it is suggested to grab it with the help of at least one other person as shown in figure

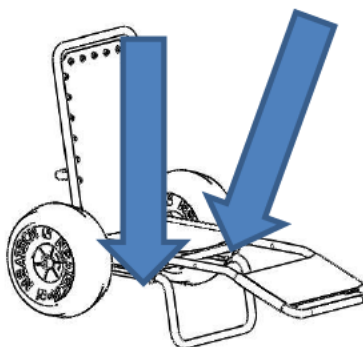


Figure 1



WARNING

It is not possible to use the product as a seat in a motor vehicle.

1 PRODUCT INFORMATION

Dear Friend,

JoB responds to the need of people with disabilities to have a specific product for the holidays and especially for bathing. Neatech has created and produced <J'amme 'o bagn> which is not only eye-catching invitation in Neapolitan to bathe in the sea, instead the JoB it has effectively solved the problem of the inability of people with disabilities or temporary limitations to enter the water, thus satisfying the natural desire of sea and atmosphere of fun and lightheartedness that characterize a vacation day. The vacation of a person with a disability should not be conceived as a simple moment of escape and entertainment, but also as an opportunity for growth.

Thanks to its waterproof wheels with variable section and controlled deformation, the JoB enables users and their companions to face pebbles, gravel, snow and dirt, transforming each new obstacle into a source of entertainment.

To facilitate the transport it has been studied a tightening system with no fixed restrictions, so the wheelchair can be assembled and dismantled in a few minutes and without the use of any tools. JoB It is the ideal solution for swimming pools


and bathing facilities not yet equipped for the stay of the disabled, as it allows you to adapt the structure to the rules in force for the removal of architectural barriers.

2 PREPARATION FOR FIRST USE

Inside the box there are: • Backrest with upholstery • Footrest with upholstery • Front support • Anti-tip bar • Main frame with upholstery and tube fasteners • N. 2 wheels • N. 10 wheel fasteners • Documents and manual Purchased accessories	 PACKAGING DISPOSAL To properly recycle the packaging materials follow instructions provided by your local waste disposal service.
---	--

Before starting to use the wheelchair please check if all described components are present. **If not, please contact as soon as possible the vendor.**

2.1 Mounting the wheelchair

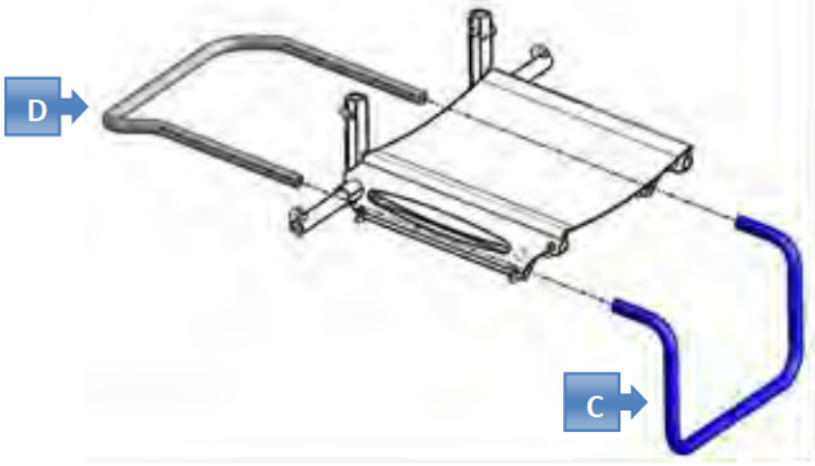	
Figure 2 • Insert anti-tip bar (D) and front support bar (C) into the main chassis (E) as shown in Figure. • Secure the two bars with fasteners. • Make sure tubes are correctly locked.	



WARNING

Tubes (D) and (C) are of the same width but (D) is longer.

If tubes (D) and (C) are inverted, the wheelchair will result unbalanced with possible risk for the user.

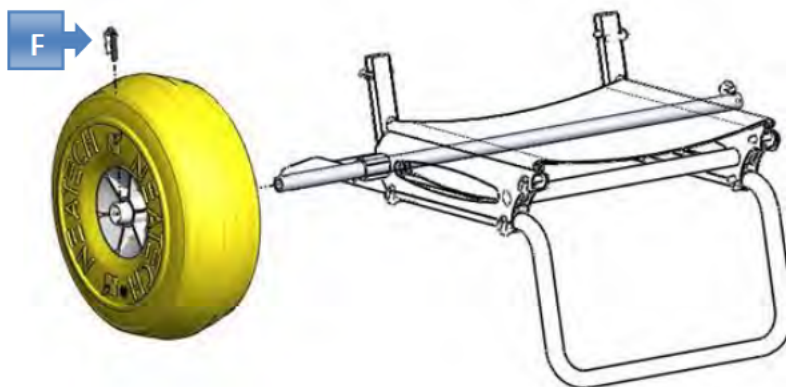


Figure 3

- Insert the wheel (lh and rh) on the axis of the chassis as shown in Figure.
- Use the special housing hole on the hub, previously aligned with the hole present on the axle of the frame and insert the fastener (F).
- Make sure each wheel fastener is correctly locked on hub.



WARNING

Make sure that the looking of the wheels and tubes are trapped in their respective housings of the frame to prevent loss.

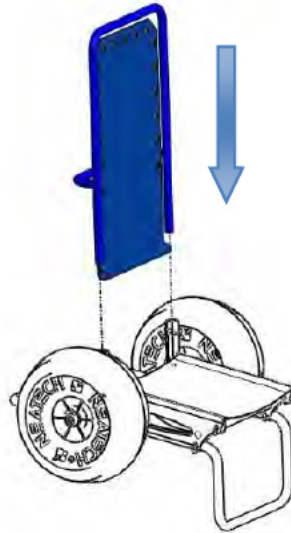


Figure 4

- Insert the backrest into the bushings on the main frame as shown in Figure.
- Insert the 2 fasteners deeply through proper holes to lock it.
- Make sure fasteners are correctly locked on main frame.



WARNING

It is very important to carefully fix the backrest with the fixing clips and the strip to avoid possible dismounting

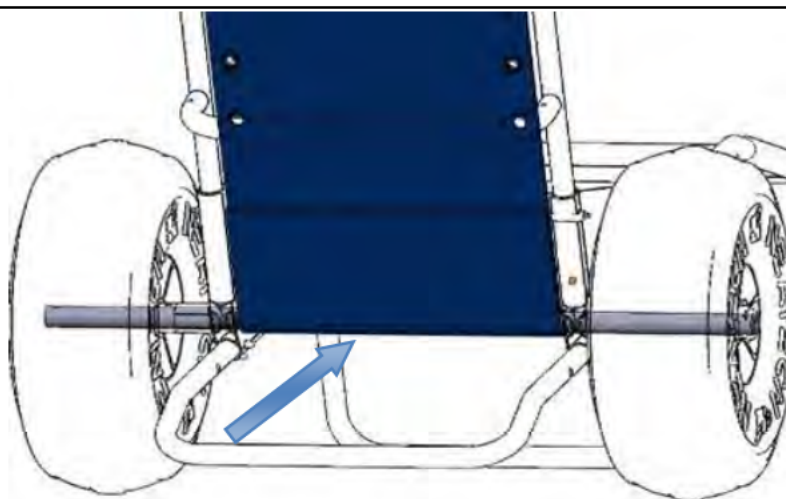


Figure 5

- Fix strip located on bottom end of the backrest, passing it around wheels' axle as shown in Figure.

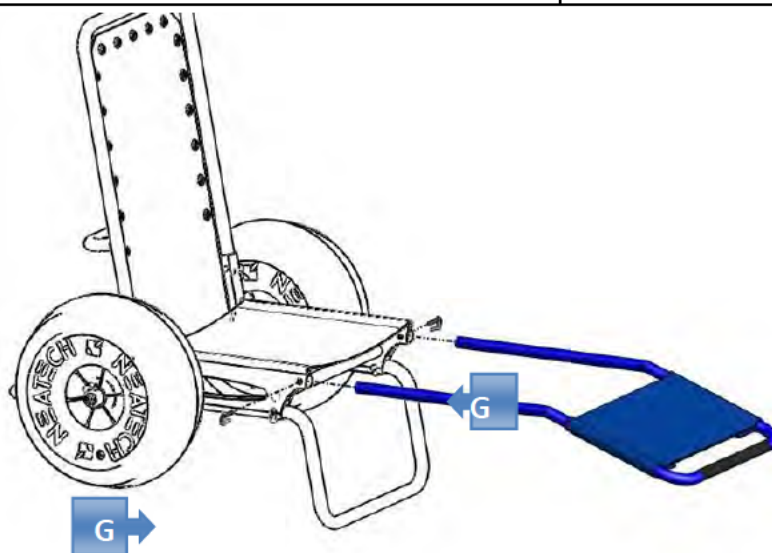


Figure 6

- Insert the legbar in the tubes of the chassis as shown in Figure and lock it with the 2 fasteners (G).



WARNING

Tighten correctly the two tubes.



Disassembly of the wheelchair

In order to disassemble the chair run the described operations in reverse order.

3 ADJUSTMENTS

Adjustment	Type of operation
Legrest adjustment	B - Assistant Operation intended to be performed by an assistant.

Table 1

3.1 Legrest adjustment



Average needed time:

<5 min

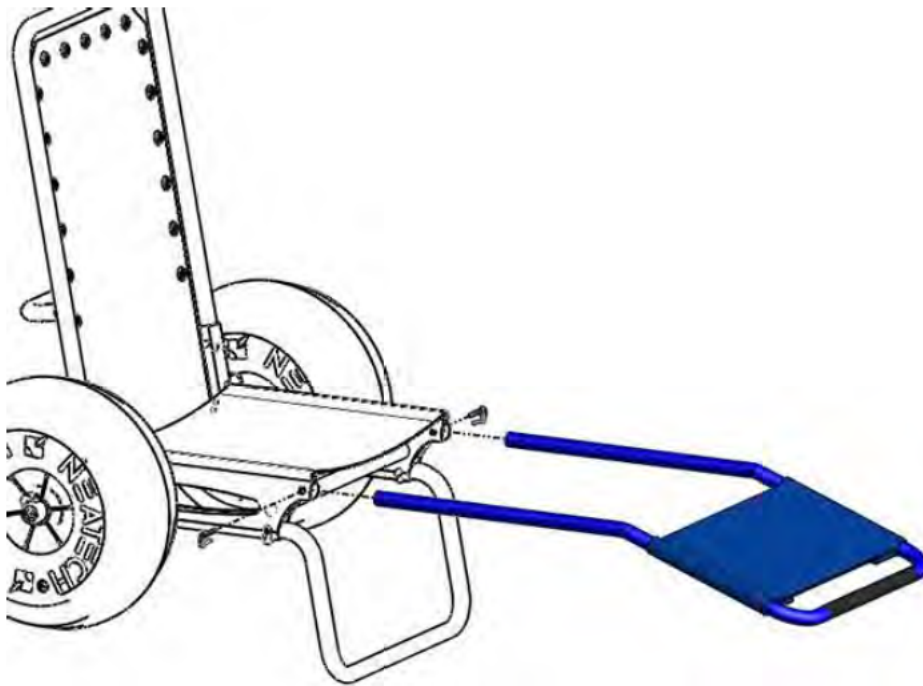


Figure 7

- Remove the fasteners as shown in Figure.
- Put the legrest in one of the three available positions.
- Insert again the fasteners.

4 USE OF THE PRODUCT

WARNING

For use in safety of the wheelchair and to avoid situations of risk for the user, verify that **the circumference of the wheels is between 1380 and 1400 mm.**

If not, please inflate or deflate the wheel to within the recommended range.



Figure 8



WARNING

Do not use wheelchair in deep or agitated water.

Presence of an assistant is obligatory when wheelchair is into the water.

4.1 Armrests

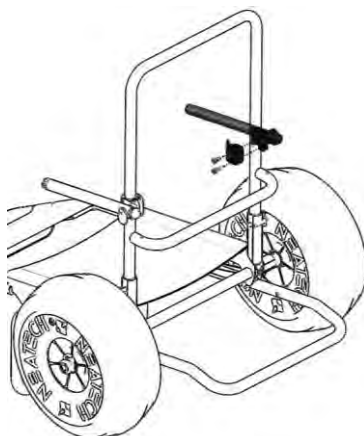


Figure 9

- Assemble the armrest (rh and lh) to the backrest as shown in Figure. To ensure a proper connection tighten the two screws shown in Figure.

ALLEN WRENCH



6 mm

4.2 Third wheel kit

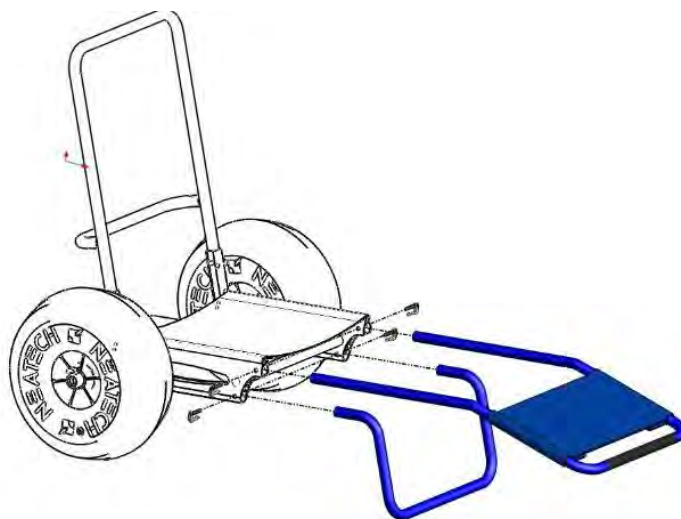


Figure 10

- Remove the leg bar and the front support bar simply acting on the fasteners as shown in Figure.

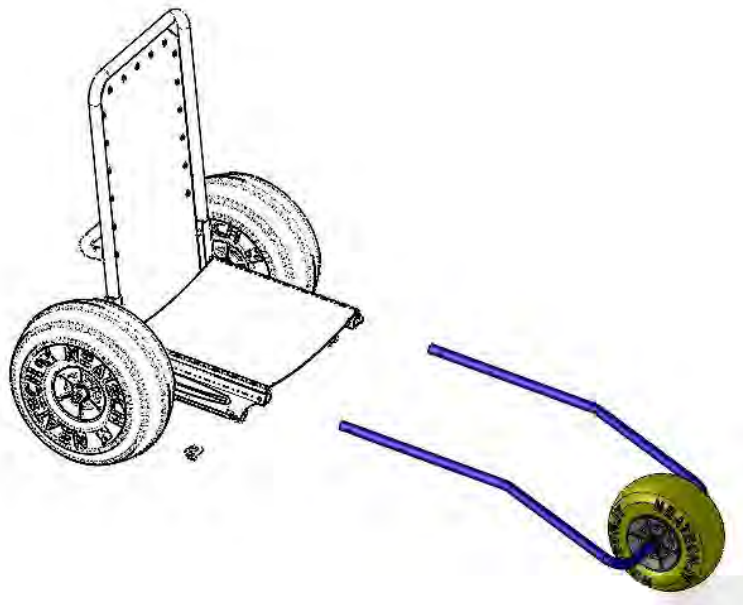


Figure 11

- Connect the third wheel kit inserting the tubes into the holes in front of the frame as shown in Figure.
- Apply the fasteners into the appropriate holes to ensure proper connection.



ADJUSTMENT

It is possible to choose among three different holes according to the desired depth.

The depth adjustment is the same adjustment as for the legbar. For more information see section 3.1.



WARNING

For use in safety of the wheelchair and to avoid situations of risk for the user, verify that **the circumference of the wheel is between 810 and 840 mm.**

If not, please inflate or deflate the wheel to within the recommended range.

4.3 Sun awning

Inside the box there are:

- Bag with awning and the straps for the coupling (1)
- Supports for the awning (2)
- Awning's support bars (3)



Figure 12



Figure 13

- From the support bars of the awnings remove the pipe cover black identifying the hole that has no outer tabs.



Figure 14

- Insert the yellow flat supports in the hole just identified and then insert the support for the awning by sliding it on the backrest. Keep the media flat-head upright and facing up.



Figure 15

- Attach the strap to the tubular on the back of the backrest so that the hook is facing towards the outside.



Figure 16

- Attach the strap on the back of the awning to the hook of the belt previously fixed to the tube and pull the belts till they are very stretched.
- Attach the ties of the awning to the tubular of the legrest.



Figure 17

- Adjust the tension of the tie rods acting on cord lock till the awning is stable with respect to the chair.

4.4 Other accessories

	Head cushion
	Carrying bag
	Inflatable collar
	Body support
	Inflatable footrest

4.5 Transport and storage

If you are not willing to use the wheelchair for a long period, keep it safe in a clean area and away from heat.

If it is necessary to transport the wheelchair to facilitate operation disassemble the wheelchair and use the bag available as accessory.

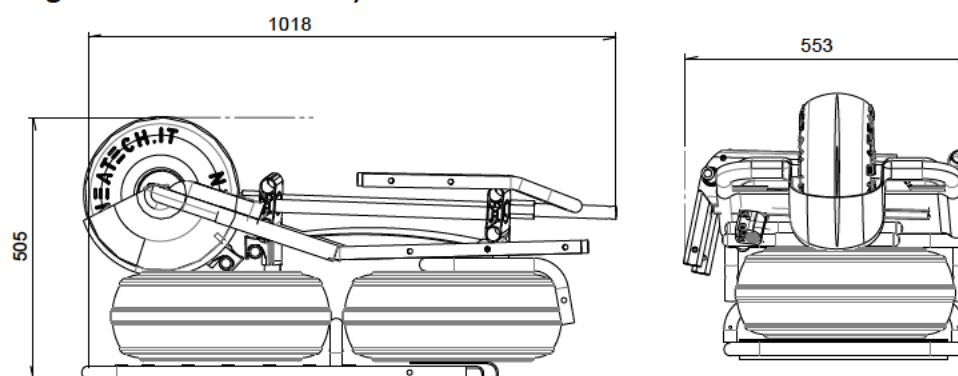


Figure 18

The wheelchair can be transported in the storage compartment of the vehicle even without any package. The unoccupied wheelchair can be also shipped via air. If it is necessary to ship the wheelchair it is very important to protect it with an appropriate package.

It is not possible to provide a universal package, so the user should provide himself for it. The used package must be water and dust resistant and strong enough to protect the wheelchair from any hurts. When inserting the wheelchair into the package protect any protruding parts with some foam or similar.

5 MAINTENANCE

Please remind that the wheelchair is intended exclusively for the carriage of seated people. Below there are some precautions for the use of the wheelchair, which it is recommend to follow, in order to ensure a safe use and a long duration.

Regular maintenance helps to keep intact the functionality and safety of the wheelchair. Inadequate or lack of care and maintenance may cause a limitation of the warranty.

- Avoid prolonged contact of the wheelchair with water. It may cause oxidation of the metal parts.
- Avoid long exposure of the wheelchair to direct sunlight.



WARNING

Any work on the wheelchair must be performed by an authorized service centre.

5.1 Cleaning

To clean the product do not use high-pressure water spray devices. For plastic and metal parts use a soft cloth dampened with mild detergent. For the upholstery, linings, seat and back covers use warm water and mild detergent. Do not use stain removers, solvents, acids, etc.

5.2 Controls to be performed on the product

Type of operation	
A	Operation intended to be performed by the user.
B	Operation intended to be performed by an assistant.
C	Operation intended to be performed by an authorized service center.

Operation	Frequency	Type of operation
Clean the wheelchair	Weekly	B - Assistant
Check if the pressure is the one indicated on tires and in section 4 - 4.2.	Weekly	B - Assistant
Check the presence of all fasteners.	Monthly	B - Assistant

5.3 Tire puncture

- | | |
|---|--|
| <ul style="list-style-type: none">• Clean damaged surface with ethylic alcohol.• When surface is dry; apply few drops of a glue for rigid PVC.• Wait some minutes to let glue dry, then inflate wheels as described in this manual. | |
|---|--|



WARNING

When finished, please be sure that there is no air leakage.

5.4 Reuse

The product is suitable for reuse. Before dispensing it, the product must be cleaned, and subjected to maintenance. The operating instructions are included in this manual and must also be provided when the product is passed on.



WARNING

This operation must be performed only at an authorized service center.

6 SPECIFICATIONS

Total mass (min-max) [kg]	~13.0 kg
Maximum user mass	125 kg
Expected lifetime of the assistive product [years]	5
Operating temperature	MIN -10 °C – MAX +40 °C
Storage temperature	MIN -20 °C – MAX +45 °C

6.1 Dimensions

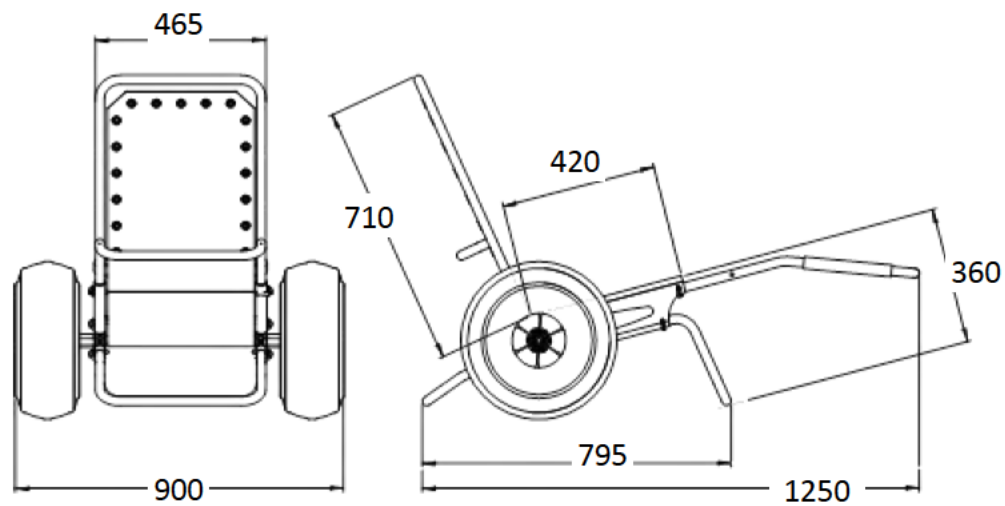


Figure 19

7 WARRANTY TERMS

Definitions

- **MANUFACTURER:** Manufacturer means the legal person who manufactures a product. For the scope of the following document, the manufacturer is:

Neatech.it s.r.l.

via Antonio de Curtis 4/A – 80040 – Cercola (NA)

Tel. +39 081 555 1946 – info@neatech.it – neatech@pec.it – www.neatech.it

P.IVA IT04812481218 – REA NA715393

- **CUSTOMER:** Customer means the natural or legal person who buys a product from the manufacturer. For the scope of the following document, the Customer is the holder of the financial document issued by the manufacturer following the supply of the product.
- **PRODUCT:** Product is the good supplied by the manufacturer to the customer in execution of a purchase order

Scope of the guarantee

The manufacturer undertakes to remedy any defect, lack of quality or lack of conformity of the products related to him as a result of design, construction errors or defects in the material that occurred during the warranty period.

Period of application of the guarantee

The warranty period begins with the customer's billing date. The duration of the warranty period varies according to the type of product.

Type of product	Warranty period [months]
Wheelchairs	24
Wheelchair parts and accessories (excluding batteries)	24

Table 2

Any repairs or replacements under warranty do not alter the original period of application of the warranty.

Exclusions

The warranty does not cover wearing parts.

Following there are some examples of parts subject to wear.

- Wheels
- Padded parts (polyurethane foams, viscoelastic)
- Upholstery (including canvas)
- Handles and knobs
- Fuses and bulbs
- Electric motor brushes
- Filters

The warranty is limited to damage related to the product and cannot under any circumstances cover damage caused to third parties as a result of product failure.

The warranty does not cover parts damaged by overload, inappropriate use, alterations and repairs made by unauthorized third parties. The warranty is not valid in the event of tampering, incorrect storage, incorrect or unauthorized maintenance.

7.1 Serial number

For any report or assistance request, please communicate the serial number mentioned on the label in the position shown in figure.

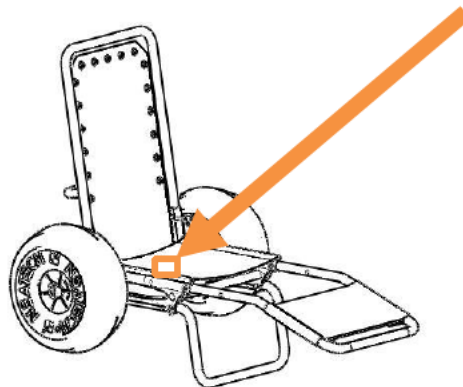


Figure 20

[illegible]

[illegible]

