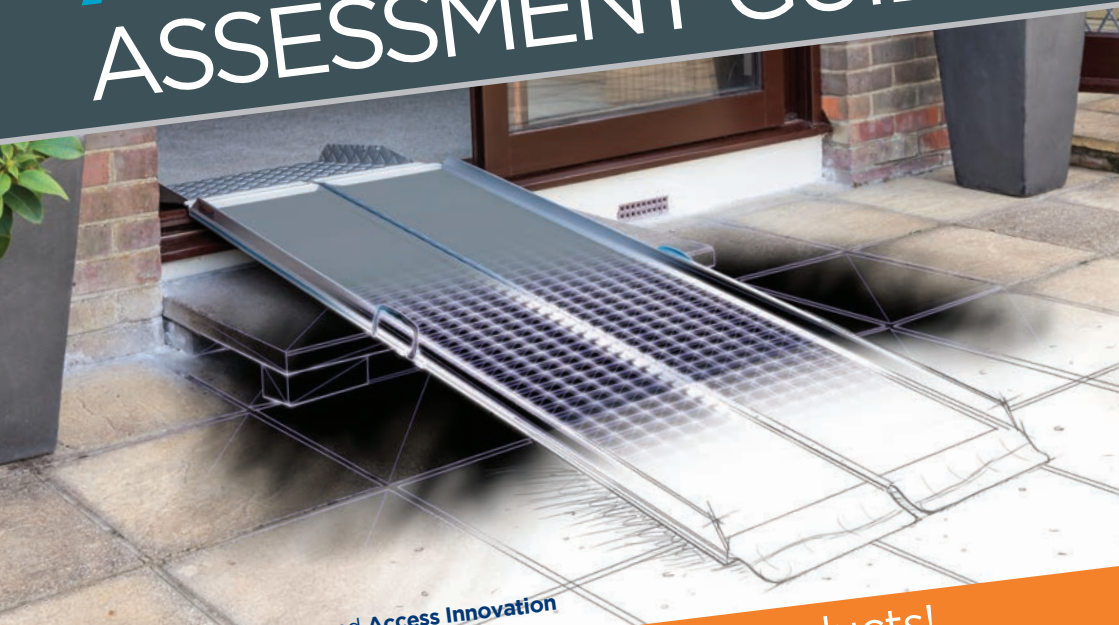


rampCENTRE®

your **Access & Ramp** ASSESSMENT GUIDE



Leading **Ramp** and **Access Innovation**

Don't Miss the **ALL NEW** Products!

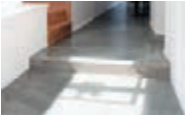



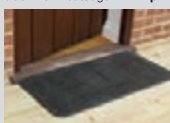


























The only Ramp Specifier Guide you'll ever need



Choosing the correct ramp is important, just like finding the right size shoe or the correct glasses. Please don't hesitate to give our experienced team a call if you need any assistance.

Choose your ramp type

 <p>Internal Steps/Thresholds</p>  <p>doorline™ Threshold Wedge p11</p>  <p>doorline™ variwedge p12</p>  <p>doorline™ neatedge p12</p>  <p>The Butterfly Ramp® p10</p> 	 <p>Steps/Stairs</p>  <p>Aerolight® p4</p>  <p>Aerolight®-High Rise p7</p>  <p>perma®-original p16</p>  <p>Welcome® p23</p> 	 <p>Wood/Level Access Door Threshold</p>  <p>Aerolight® p4</p>  <p>doorline™ p10</p>  <p>perma®-original p16</p>  <p>Welcome® p23</p> 	 <p>uPVC Door Threshold</p>  <p>Aerolight®-Up&Over p7</p>  <p>The Butterfly Ramp® p10</p>  <p>perma®-adjust p17</p>  <p>Welcome® p23</p> 	 <p>Vehicle Access</p>  <p>Aerolight®-Lifestyle p5</p>  <p>Aerolight®-Broadfold p6</p>  <p>ultra®-telescopic p14</p> 
--	--	--	---	--

Key icons



Suitable for
all weathers



Compact
Storage



Lightweight
Solution



Look out for
great tips



Leave
in-place



Fast
install

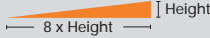
How steep?

1:6 (10°)*



Acceptable in some cases for short threshold ramps

1:8 (8°)*




A useful balance between gradient and space suitable for portable ramps in some assisted applications and events, etc.

1:12 (5°)



Often required by social services/local authorities; minimum required by Building Regulations Document M for residential installed ramps

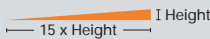1:15 (4°)



Sometimes required by social services/local authorities for independent access, but typically for installed public access ramps

*Assess safety/risk

How long?

Try the
ramp gradient
calculator on
our website!

- a. Measure the step height/level difference.
(Total height of steps or height to resting point on vehicle)



Remember:
"The longer the ramp,
the easier the slope"




- b. Use the **Ramp Length Guide**

Height x Gradient = Length Required

Height (up to)	Length		
	1:6	1:8	1:12
10cm (4ins)	60cm (2ft)	80cm (2ft 8ins)	120cm (4ft)
15cm (6ins)	90cm (3ft)	120cm (4ft)	180cm (6ft)
20cm (8ins)	120cm (4ft)	160cm (5ft 4ins)	240cm (8ft)
25cm (10ins)	150cm (5ft)	200cm (6ft 8ins)	300cm (10ft)
30cm (1ft)	180cm (6ft)	240cm (8ft)	360cm (12ft)
45cm (1ft 6ins)	270cm (9ft)	360cm (12ft)	540cm (18ft)
60cm (2ft)	360cm (12ft)	480cm (16ft)	720cm (24ft)

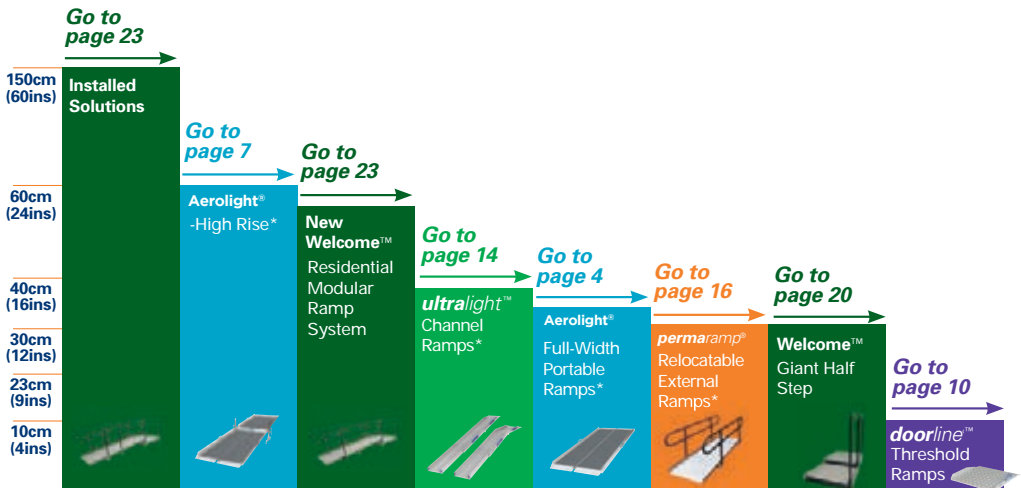Select by height range



*Heights shown are based on 1:8 gradient.

rampCENTRE® **NEW Aerolight®** Full-width portable ramps



So **VERSATILE**, So **EASY-TO-USE**,
 So **QUICK TO DEPLOY**, So **SAFE**,
 So **STRONG**, Now Even **LIGHTER!**

- Great for all types of wheelchairs, rollators and scooters (2, 3, 4, 5 and 6 wheel types)
- Easy to see and walk on – with **NEW** reflective lip feature and full-width grip surface
- Available with different methods of folding to suit your needs
- **NEW** unique tapered lip design to fit narrow doorways

NEW Aerolight®-Xtra Premium folding suitcase ramp



A proven and popular workhorse



Custom sizes available



7 sizes

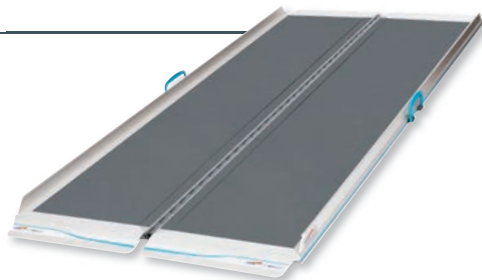


2ft - 8ft long



Up to 450kg capacity

- ▶ Folds and clips for convenient, compact storage
- ▶ Flat underneath for the best possible step clearance
- ▶ Very easy to carry and deploy
- ▶ Anti-slip surface folds inside when stored and carried for protection and safety
- ▶ Contrasting colour edge strip increases visibility



Code	Length	Capacity	Weight
AX6	60cm (2ft)	450kg (71st)	4.9kg (11lbs)
AX9	90cm (3ft)	450kg (71st)	6.6kg (15lbs)
AX12	120cm (4ft)	400kg (63st)	8.6kg (19lbs)
AX15	150cm (5ft)	350kg (55st)	10.4kg (23lbs)
AX18	180cm (6ft)	350kg (55st)	12.3kg (27lbs)
AX21	210cm (7ft)	300kg (47st)	14.2kg (31lbs)
AX24	240cm (8ft)	250kg (39st)	16kg (35lbs)

Ramp width (excluding handles) 76cm (30ins). Lip width 71.4cm (28.1ins).

NEW Aerolight®-Classic

Premium non-folding ramp

So easy – simply place and use



- Great portable solution for thresholds, kerbs and steps
- Well suited for public access
- Perfect for walking users


Code	Length	Capacity	Weight
AC6	60cm (2ft)	400kg (63st)	4.3kg (9lbs)
AC9	90cm (3ft)	350kg (55st)	5.6kg (12lbs)
AC12	120cm (4ft)	350kg (55st)	7.1kg (16lbs)
AC15	150cm (5ft)	350kg (55st)	8.5kg (19lbs)

Ramp width (excluding handles) 76cm (30ins). Lip width 71.4cm (28.1ins).

NEW Aerolight®-Lifestyle

Premium multi-folding ramp

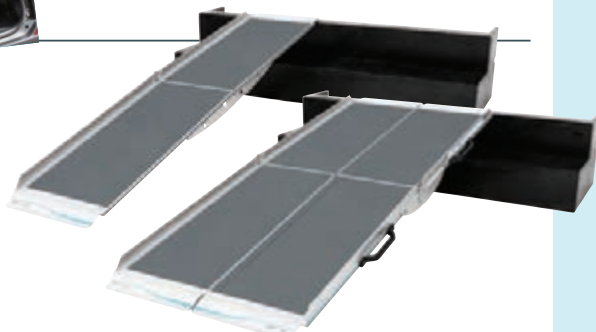Full-width safety and easy handling; great for out-and-about access and home use!



- Safety hinges with pinch protection
- Highest strength-to-weight ratio scooter ramp available
- Easy vehicle loading
- Two halves make it very light to carry and easy to store

Code	Length	Capacity	Weight (Each Half)	Stored Dimensions L x W x D
AL18	180cm (6ft)	350kg (55st)	7.8kg (17lbs)	95 x 37 x 25cm (37 x 15 x 10ins)
AL21	210cm (7ft)	300kg (47st)	8.8kg (19lbs)	110 x 37 x 25cm (43 x 15 x 10ins)
AL24	240cm (8ft)	250kg (39st)	9.8kg (22lbs)	125 x 37 x 25cm (49 x 15 x 10ins)
AL27	270cm (9ft)	250kg (39st)	10.9kg (24lbs)	140 x 37 x 25cm (55 x 15 x 10ins)

Ramp width (excluding handles) 74cm (29ins). Lip width 69.6cm (27.4ins).

NEW Aerolight® - Full-width portable ramps

NEW Aerolight®-Broadfold

Premium folding ramp

Quick, easy & safe;
great for vehicle
access



A simple and effective scooter ramp, the NEW Aerolight®-Broadfold folds to half-length for compact storage.

- Now with 4 rigid handles for safety and ease-of-use
- Easy vehicle loading - superb solution for taxi and minibus ramp applications
- Highest capacity broadfold ramp on the market!
- So easy to deploy!
- Full-width lip locates firmly into vehicle




Code	Length	Capacity	Weight	Stored Dimensions L x W x D
AB15	150cm (5ft)	350kg (55st)	12.1kg (27lbs)	80 x 19 x 76cm (31 x 7 x 30ins)
AB18	180cm (6ft)	350kg (55st)	13.5kg (30lbs)	95 x 19 x 76cm (37 x 7 x 30ins)
AB21	210cm (7ft)	300kg (47st)	14.9kg (33lbs)	110 x 19 x 76cm (43 x 7 x 30ins)
AB24	240cm (8ft)	250kg (39st)	16.4kg (36lbs)	125 x 19 x 76cm (49 x 7 x 30ins)

Ramp width (excluding handles) 76cm (30ins). Lip width 71.4cm (28.1ins).

Aerolight®-Travel Kerb Ramp



Extremely compact and lightweight

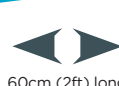A simple ready-to-go
solution



The all-new Aerolight-Travel Kerb Ramp is your quick and easy to use answer to out-and-about access.

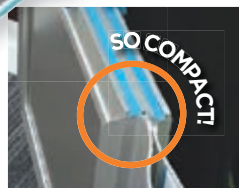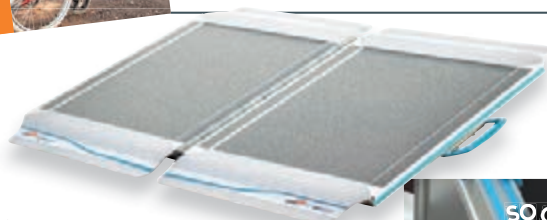- Just right for kerbs and steps
- Handy carry bag to fit on the back of the wheelchair
- Great for any mobility device
- Full-width anti-slip surface for safety



Code	Length	Ramp Width	Lip Width	Weight	Capacity
AXT	60cm (2ft)	77cm (30ins)	71.4cm (28.1ins)	4.5kg (10lbs) (excluding bag)	350kg (55st)

Kerb Height 7cm (3ins)

NEW Aerolight®-High Rise

NEW Aerolight® combinations up to 16ft long

Achieve incredible lengths and heights with lightweight, portable ramps

- Fully adjustable for any height and gradient
- Ramp joining support (AJ4) is also available separately
- Save £'s by buying as a kit



NEW Aerolight®-Xtra

AJ4 joining support

NEW Aerolight®-Xtra

NEW Aerolight®-High Rise kit

Code	Includes	Length	Capacity
HR30	1 x AJ4, 2 x AX15	300cm (10ft)	350kg (55st)
HR36	1 x AJ4, 2 x AX18	360cm (12ft)	350kg (55st)
HR42	1 x AJ4, 2 x AX21	420cm (14ft)	300kg (47st)
HR48	1 x AJ4, 2 x AX24	480cm (16ft)	250kg (39st)

Code	Description	Height	Capacity	Weight
AJ4	Joining Support: Includes easy-use handle fixings for height adjustment. Can be used with any NEW Aerolight® ramps.	11-45cm (4-18ins)	350kg (55st)	4.7kg (10lbs)

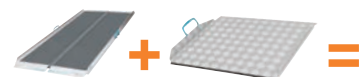
Ramp width (excluding handles) 76cm (30ins). Lip width 71.4cm (28.1ins).

NEW Aerolight®-Up&Over

NEW Aerolight®-Xtra and threshold ramp combinations

The complete up-and-over solution for thresholds - great for uPVC thresholds

- Versatile, quick and easy to deploy
- Solves most threshold types and heights
- Portable and easy to store
- Seamless transition right over the threshold
- Save £'s by buying as a kit



NEW Aerolight®-Xtra

Doorline-Multi (page 10)

NEW Aerolight®-Up&Over kit

Code	Includes	External Ramp Length	Capacity
UO12	1 x DM5, 1 x AX12	120cm (4ft)	300kg (47st)
UO15	1 x DM5, 1 x AX15	150cm (5ft)	300kg (47st)
UO18	1 x DM5, 1 x AX18	180cm (6ft)	300kg (47st)
UO21	1 x DM5, 1 x AX21	210cm (7ft)	300kg (47st)
UO24	1 x DM5, 1 x AX24	240cm (8ft)	250kg (39st)

NEW Aerolight®-Accessories



Aerolight® Wall Strap Kit Code: AWK

Quick-fit strap holds any Aerolight® ramp vertically against wall for tidy storage.



Aerolight® Joining Support Code: AJ4

Use to join 2 ramps together.



Aerolight®-Connect Platform

For use with Aerolight ramps – excellent for events, up to 38cm (15ins) high

CP12

Aerolight-Connect platform 122 x 122cm c/w safety kerbs, doorplate, adjustable legs and rubber feet (no handrails)

CPHK

WRSK

WRSK1

WRSK2

Single handrail kit for above Step kit for above (no handrails)

Step handrail kit for one step

Step handrail kit for two steps



Aerolight® Upright Storage Bracket Kit Code: ASBV

Compact vertical storage for Aerolight-Xtra ramps. Comes with locking top section for extra security. (Max. 180cm ramp).



Aerolight® Horizontal Storage Bracket Kit - Code: ASB

Easy, ergonomic storage for Aerolight®-Classic and Xtra ramps (with padlock security).



Aerolight® Ramp Locating Peg Kit Code: ARLS

Can be fitted to any Aerolight® ramp to securely locate on upper step or in vehicle, etc. Kit includes two studs, drill bit and fixings. Underside of lip shown with locating studs.

Utilityramp™

Industrial/Commercial Portable Ramp

Rugged, yet easy to use
and portable



12 sizes



Up to 240cm-8ft long

- Great for wheelchair access or industrial use
- Unique Gripdeck™ highly anti-slip surface retains grip even when wet
- Up to 400kg weight capacity
- Standard 76cm wide or Superwide 90cm wide versions available

Code	Length	Width	Capacity	Weight
UR9*	90cm (3ft)	76cm (30ins)	400kg (63st)	9kg (20lbs)
UR12*	120cm (4ft)	76cm (30ins)	400kg (63st)	12.5kg (28lbs)
UR15	150cm (5ft)	76cm (30ins)	400kg (63st)	15.5kg (34lbs)
UR18	180cm (6ft)	76cm (30ins)	400kg (63st)	19kg (42lbs)
UR21	210cm (7ft)	76cm (30ins)	350kg (55st)	22.5kg (50lbs)
UR24	240cm (8ft)	76cm (30ins)	300kg (47st)	25.5kg (56lbs)
URW9*	90cm (3ft)	90cm (35ins)	400kg (63st)	10.5kg (23lbs)
URW12*	120cm (4ft)	90cm (35ins)	400kg (63st)	14.5kg (32lbs)
URW15	150cm (5ft)	90cm (35ins)	400kg (63st)	18.5kg (41lbs)
URW18	180cm (6ft)	90cm (35ins)	400kg (63st)	22kg (49lbs)
URW21	210cm (7ft)	90cm (35ins)	350kg (55st)	26kg (57lbs)
URW24	240cm (8ft)	90cm (35ins)	300kg (47st)	30kg (66lbs)

* No wheels. Kerb Height 7cm (3ins)



Rollout-Trackway™ & Rollout-Trackway Patio™

Instant access solution for grass, gravel, and other uneven surfaces



Rollout-Trackway – Pathway Solutions

- Great for events and short-term uses, and can be left in place for long-term use
- Totally weather-resistant with full drainage - perfect for the unpredictable British weather
- Instant deployment for immediate use
- Extremely durable, and can cover grass, gravel and cobbles etc
- Saves expensive and invasive construction work
- Multiple pieces can be used for larger areas

Rollout-Trackway™

Code	Length
RT15	150cm (4ft 11ins)
RT25	250cm (8ft 2ins)
RT50	500cm (16ft 5ins)

Width 90cm (35.4ins)

Rollout-Trackway Patio™

Code	Length
RTP20	2 x 1.8m (6ft 7ins x 5ft 10ins)

rampCENTRE[®] doorline[™] Threshold ramps

Whatever the threshold, the Doorline range can cope.

Huge range • **Adjustable options**
Easy-to-use • **Excellent finish**

- Very simple to deploy
- Can be used either or both sides of a door
- Removable and leave-in-place solutions available

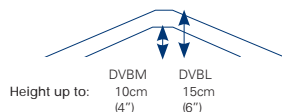


The *ButterflyRamp*[®]

This totally unique bridge ramp automatically adjusts as it is placed to suit all threshold types, with no user adjustment required!

Your benefits:

- So fast and easy to use, and perfect for any mobility device
- Low profile design for a shallow and safe gradient
- Two convenient sizes to cover thresholds of up to 6" high
- All-weather full-width anti-slip surface folds inside when carrying for extra safety, with integral comfort carrying handles
- Registered designs, patents and patent pending



Code	Height Range From	Height Range To	Width	Overall Length	Capacity	Weight	Stored Dimensions L x H x D
DVBM	Zero	10cm plus (4ins plus)	70cm (28ins)	90cm (36ins)	350kg (55st)	6.8kg (15lbs)	70 x 48 x 4cm (28 x 19 x 1.5ins)
DVBL	Zero	15cm plus (6ins plus)	70cm (28ins)	120cm (48ins)	300kg (47st)	9.1kg (20lbs)	70 x 63 x 4cm (28 x 25 x 1.5ins)

doorline[™]-multi

Your benefits:

- Anodised finish for many years of service
- Comfortable, foldaway, rattle-free handle
- Designed for use with our NEW Aerolight and Ultralight ramps to easily bridge any threshold
- Unique lip design with rubber anti-slip to safely eliminate threshold damage
- Extremely lightweight for quick and easy use, with full 300kg capacity



Code	Max Recom. Height	Total Length	Capacity	Weight
DM5	10cm (3.9ins)	56cm (22ins)	300kg (47st)	3.6kg (7.8lbs)
DM7	13cm (5.1ins)	74cm (29ins)	300kg (47st)	4.7kg (10.4lbs)

Width 70cm (28ins)

“Designed to overcome
uPVC door thresholds”

doorline™-Threshold Wedge



“Simple indoor and outdoor threshold solution”



Your benefits:

- 6 popular inch sizes immediately available from stock
- Superb alternative to fibreglass wedge ramps
- Fully sanitisable and washable with pressure and heat
- Unique new robust and lighter weight design
- Excellent all-weather grip surface



Easily fixed in place for extra security

Code	Replaces Old Code	Height	Length	Width	Capacity	Weight
DTW1.5	DW4	4cm (1.5ins)	20cm (7.8ins)	74cm (29ins)	350kg (55st)	1.7kg (3.7lbs)
DTW2	DW6	5cm (2ins)	27cm (10.6ins)	74cm (29ins)	350kg (55st)	2.3kg (5lbs)
DTW3	DW8	7.5cm (3ins)	40cm (15.7ins)	74cm (29ins)	350kg (55st)	4.5kg (10lbs)
DTW4	DW10	10cm (4ins)	50cm (19.7ins)	74cm (29ins)	350kg (55st)	5.5kg (12lbs)
DTW5	-	12.5cm (5ins)	70cm (27.6ins)	74cm (29ins)	350kg (55st)	7.3kg (16lbs)
DTW6	-	15cm (6ins)	90cm (35.4ins)	74cm (29ins)	350kg (55st)	9.2kg (20.2lbs)

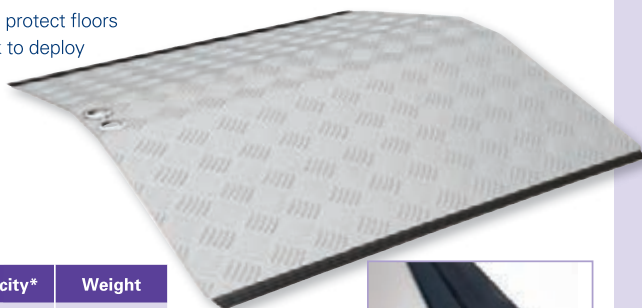
doorline™-bridge



Your benefits:

- Unique comfortable carrying eyelets for safety and comfort
- High quality finish to preventing stains and corrosion for years of use
- Rubber strips to prevent slipping and protect floors
- Extremely lightweight and very quick to deploy

“Bridge almost any threshold!”



Code	Height	Length	Capacity*	Weight
DB4	4cm (1.6ins)	50cm (20ins)	300kg (47st)	3kg (6.6lbs)
DB6	6cm (2.4ins)	58cm (23ins)	300kg (47st)	3.5kg (7.7lbs)
DB8	8cm (3.1ins)	75cm (29ins)	300kg (47st)	5.4kg (11.9lbs)
DB10	10cm (3.9ins)	100cm (39ins)	300kg (47st)	7.6kg (16.8lbs)

Width 70cm (28ins) *Capacity assumes ramp rests on threshold



Gripline™ strips to prevent slipping and protect floor

doorline™ - versatile threshold ramps

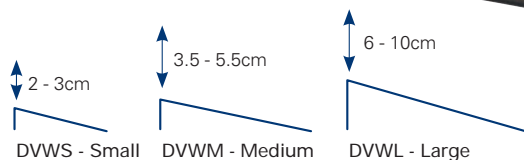
doorline™-variowedge



Adjust height for any threshold up to 10cm (4ins)


Your benefits:

- Extremely durable designer-finish to suit any internal décor
- Match threshold height to within the thickness of a coin (adjusts in 5mm increments)
- Can be fixed in place for maximum stability and security
- Rubber strip to grip and protect floor



Code	Length	Weight
DVWS	15cm (5.9ins)	1.3kg (2.9lbs)
DVWM	25cm (9.8ins)	2.1kg (4.6lbs)
DVWL	50cm (19.7ins)	6.1kg (13.5lbs)

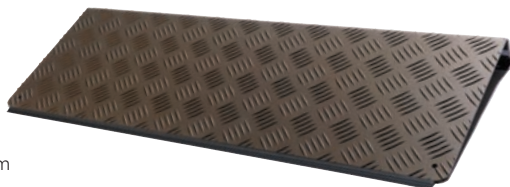Width 74cm (29.1ins). Capacity 300kg (47st).




Assemble or adjust for re-use in less than a minute!

doorline™-neatedge



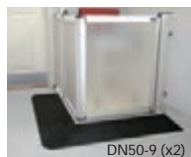
Unobtrusive, strong and versatile wedge ramps up to 152cm (5ft) wide



DN50-9 (x2)



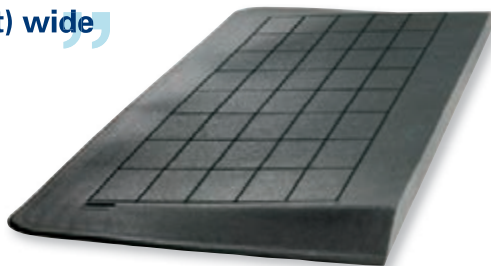Chamfered edges (not DN6) allow easy access from any angle



Code	Height	Length	Width	Weight
DN15-9	1.5cm (0.6ins)	12cm (4.5ins)	90cm (35ins)	1.4kg (3lbs)
DN25-9	2.5cm (1ins)	20cm (8ins)	90cm (35ins)	3.1kg (6.8lbs)
DN35-9	3.5cm (1.4ins)	28cm (11ins)	90cm (35ins)	5.9kg (13lbs)
DN50-9	5cm (2ins)	40cm (16ins)	90cm (35ins)	11kg (24.3lbs)
DN75-9	7.5cm (3ins)	60cm (24ins)	90cm (35ins)	18kg (39.7lbs)
DN100-9	10cm (4ins)	80cm (31.5ins)	90cm (35ins)	31kg (68.3lbs)
DN125-9	12cm (5ins)	100cm (39.4ins)	90cm (35ins)	51kg (112.4lbs)
DN2	2.2cm (0.9ins)	20cm (8ins)	107cm (42ins)	3.8kg (8.4lbs)
DN3	3.2cm (1.3ins)	31cm (12ins)	107cm (42ins)	7.3kg (16.1lbs)
DN4	4.4cm (1.8ins)	46cm (18ins)	127cm (50ins)	17kg (37.5lbs)
DN5	5.1cm (2.0ins)	53cm (21ins)	152cm (60ins)	28kg (61.7lbs)
DN6*	5.7cm (2.3ins)	75cm (30ins)	76cm (29.5ins)	21.5kg (47.4lbs)
DNF	Adhesive for permanent fixing (290ml), including application gun			

*DN6 does not have chamfered edges. Can be trimmed or cut to size as required. Capacity 500kg (79st)

doorline™-neatslope



doorline™ - versatile threshold ramps

Eliminate low threshold trip hazards!

Your benefits:

- Available in 3 useful heights
- Easy to trim if required
- Available with or without self-adhesive backing




Code	Height	Length	Self Adhesive Y/N
DNS6	0.6cm (0.25ins)	6cm (2.4ins)	N
DNS6A	0.6cm (0.25ins)	6cm (2.4ins)	Y
DNS10	1cm (0.4ins)	10cm (3.9ins)	N
DNS10A	1cm (0.4ins)	10cm (3.9ins)	Y
DNS16	1.6cm (0.6ins)	15cm (5.9ins)	N
DNS16A	1.6cm (0.6ins)	15cm (5.9ins)	Y

Width 90cm (35.4ins). Capacity 300kg (47st).



doorline™-custom

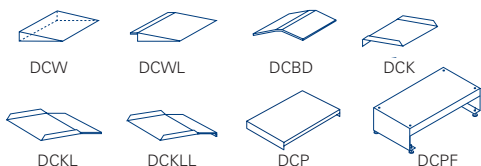Specially made to your specifications



DCW

DCBD



These ramps are individually manufactured to your specified sizes.

Choose from 8 styles for the very few applications our stock products cannot solve.

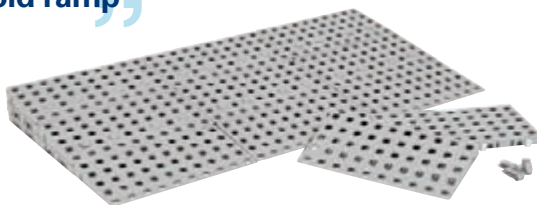
Build-a-Slope™



Totally modular threshold ramp

Your benefits:

- Easily configured for any application
- Fully weatherproof
- Quick to assemble and install
- Choose from 4 popular kits
- Compatible with other brands



Code	Max. Height	Length
BSK1	5.5cm (2.2ins)	25cm (10ins)
BSK2	9.1cm (3.6ins)	50cm (20ins)
BSK3	12.7cm (5ins)	75cm (29.5ins)
BSK4	16.3cm (6.4ins)	100cm (39ins)

Width 75cm (29.5ins). Capacity 450kg (71st).



rampCENTRE® ultralight™



Exceptionally well crafted in anodised aluminium, the Ultralight original perforated channel ramps offer an outstanding combination of compact storage and ease of use.

Anodised for durability • Unique built-in grip surface • Foldaway handles • Wide channels

- Great for narrow doorways
- Commonly used for residential and vehicle access
- Great for standard wheelchairs and four-wheel scooters

- Very compact for storage
- Extremely portable
- Great for cross-slopes


ultralight™ -folding

Ideal for those on the move



Code (Pair)	Length	Capacity (Pair)	Weight (Each Half)
PF11	110cm (3ft 7ins)	300kg (47st)	3.6kg (8lbs)
PF15	150cm (4ft 11ins)	300kg (47st)	4.6kg (11lbs)
PF20	200cm (6ft 7ins)	300kg (47st)	6kg (14lbs)

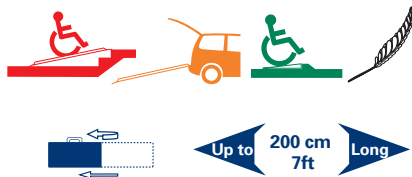Internal channel width 21cm (8ins)



ultralight™ -telescopic

Unlimited adjustability



Code (Pair)	Length		Width (Internal Channel)	Capacity (Pair)	Weight (Each Half)
	Max	Min			
PT12-2	120cm (3ft 11ins)	75cm (2ft 5ins)	19cm (8ins)	400kg (63st)	3.6kg (8lbs)
PT20-2	202cm (6ft 8ins)	117cm (3ft 10ins)	19cm (8ins)	300kg (47st)	5.8kg (13lbs)
PT20-3	199cm (6ft 6ins)	87cm (2ft 10ins)	18cm (7ins)	300kg (47st)	6.1kg (14lbs)
PT29-3	288cm (9ft 5ins)	117cm (3ft 10ins)	18cm (7ins)	250kg (39st)	8.2kg (18lbs)

3
Section

ultralight™-combi

“Amazing versatility for straight or up-and-over use”

Your benefits:

- Low kerbs and long bevel sides greatly increase clearance for scooters and power chairs
- Good choice for bridging a threshold with unequal levels each side



Code (Pair)	Length		Width (Internal Channel)	Capacity (Pair)	Weight (Each Half)
	Max	Hinged Section			
PC21-3	207cm (6ft 9ins)	75cm (30ins)	19cm (8ins)	300kg (47st)	6.5kg (14lbs)
PC28-3	276cm (9ft 1ins)	100cm (39ins)	19cm (8ins)	250kg (39st)	8.4kg (19lbs)
PCW21-3	205cm (6ft 8ins)	75cm (30ins)	25cm (10ins)	300kg (47st)	8.2kg (18lbs)
PCW28-3	294cm (9ft 8ins)	106cm (42ins)	25cm (10ins)	250kg (39st)	10.9kg (24lbs)

Extra Width

Don't miss the NEW PRODUCTS!



Aerolight-Travel Ramp
See Page 6



Butterfly Ramp
See Page 10



permaramp™-quickrails™
See Page 18

ultralight™ - compact lightweight channel ramps

ramPCENTRE® permaramp®

Relocatable external ramps

A concept that totally transforms wheelchair access, the Permaramp range is designed to be easily adjusted and fixed down or left in place for every type of threshold.

Cut-away kerbs • Easy installation • Extremely durable

- Great for semi-permanent access
- Avoid big installation projects and construction mess
- Can be left in place as long as required then re-used elsewhere
- Quick and strong access solutions
- Unique Gripdeck™ surface for incredible grip and sustained outside use



Gripdeck™
technology

Incredible Grip



permaramp®-original

Your benefits:

- Excellent for wooden door thresholds
- Contrasting colour safety edge strips increase visibility
- Superwide 90cm width version available as well as standard 76cm width: perfect for single outward-opening doors, and a great bariatric ramp solution
- Fixing options at top and bottom for secure location
- Well suited for scooters, mobility devices and sports wheelchairs
- Up to 400kg (63 stone) capacity!



“Instant installed ramp solution available in 2 widths and 6 lengths”



PO15

EXTRA
WIDTH


Code	Length	Width	Capacity
PO9	90cm (3ft)	76cm (30ins)	400kg (63st)
PO12	120cm (4ft)	76cm (30ins)	400kg (63st)
PO15	150cm (5ft)	76cm (30ins)	400kg (63st)
PO18	180cm (6ft)	76cm (30ins)	400kg (63st)
PO21	210cm (7ft)	76cm (30ins)	350kg (55st)
PO24	240cm (8ft)	76cm (30ins)	300kg (47st)
POW9	90cm (3ft)	90cm (35.4ins)	400kg (63st)
POW12	120cm (4ft)	90cm (35.4ins)	400kg (63st)
POW15	150cm (5ft)	90cm (35.4ins)	400kg (63st)
POW18	180cm (6ft)	90cm (35.4ins)	400kg (63st)
POW21	210cm (7ft)	90cm (35.4ins)	350kg (55st)
POW24	240cm (8ft)	90cm (35.4ins)	300kg (47st)

permaramp®-adjust

Your benefits:

- Easy height adjustment and installation from above
- Designed to clear protruding thresholds and sills
- Fix to any firm flat or sloping surface
- Telescopic legs can be adjusted to accommodate different threshold types, and give total height adjustment from 7.5cm (3ins) to 25cm (10ins) high
- Entry flap accessory available to transition right into the property (see below)
- All lengths available in Superwide 90cm option: a great bariatric ramp solution and perfect for single outward-opening doors, and extra-wide mobility devices

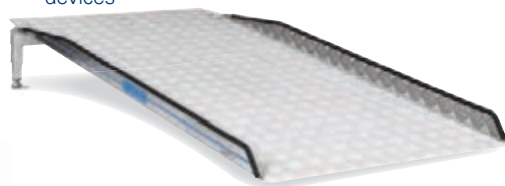Totally adjustable and self-supporting, great for uPVC thresholds



Code	Length	Width	Capacity
PA9	90cm (3ft)	76cm (30ins)	400kg (63st)
PA12	120cm (4ft)	76cm (30ins)	400kg (63st)
PA15	150cm (5ft)	76cm (30ins)	400kg (63st)
PA18	180cm (6ft)	76cm (30ins)	400kg (63st)
PA21	210cm (7ft)	76cm (30ins)	350kg (55st)
PA24	240cm (8ft)	76cm (30ins)	300kg (47st)
PAW9	90cm (3ft)	90cm (35.4ins)	400kg (63st)
PAW12	120cm (4ft)	90cm (35.4ins)	400kg (63st)
PAW15	150cm (5ft)	90cm (35.4ins)	400kg (63st)
PAW18	180cm (6ft)	90cm (35.4ins)	400kg (63st)
PAW21	210cm (7ft)	90cm (35.4ins)	350kg (55st)
PAW24	240cm (8ft)	90cm (35.4ins)	300kg (47st)

EXTRA WIDTH

permaramp®-entry flap

The ultimate semi-permanent uPVC threshold bridging accessory

Your benefits:

- Flap accessory for Permaramp-Adjust (see above)
- Hand strap for fast and easy operation
- Excellent for sliding and inward-opening doors
- Unique Gripdeck™ surface on both sides



Code	Length	Width	Capacity
PE5	45cm (18ins)	72cm (28.3ins)	300kg (47st)

NOTE: Only for use with standard width Permaramp-Adjust, ordered separately

rampCENTRE® permaramp®-quickrails™

The all-new permaramp®-quickrails™ are your answer for an instant, fixed ramp solution.

This unique range is designed to be extremely universal and easy to fit.



NEW

Pre-assembled handrails

- One side, both sides (or no sides!)
- Interchangeable left or right
- Retro-fittable
- 4 colour options



Choose from

permaramp®-adjust
totally adjustable and
self-supporting, great
for uPVC thresholds.

or

permaramp®-original
for fixing straight
to the threshold, ideal for
wooden door thresholds

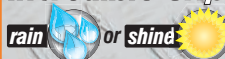
Factory-fitted handrail sockets

Super-Fast and Convenient Assembly!

No-Tools assembly
kit option

Gripdeck™
technology

Incredible Grip



2 widths

- 76cm (30ins)
- 90cm (36ins) Superwide

5 lengths

- 120cm (4ft)
- 150cm (5ft)
- 180cm (6ft)
- 210cm (7ft)
- 240cm (8ft)



Seen the
video?





**permaramp®-original
with quickrails™**

**permaramp®-adjust
with quickrails™**

Length	 One Side (FIT LEFT OR RIGHT)	 Both Sides	 One Side (FIT LEFT OR RIGHT)	 Both Sides
--------	---	--	---	--

76cm (30ins) width





120cm (4ft)	PO12R1	PO12R2	PA12R1	PA12R2
150cm (5ft)	PO15R1	PO15R2	PA15R1	PA15R2
180cm (6ft)	PO18R1	PO18R2	PA18R1	PA18R2
210cm (7ft)	PO21R1	PO21R2	PA21R1	PA21R2
240cm (8ft)	PO24R1	PO24R2	PA24R1	PA24R2

90cm (36ins) Superwide

120cm (4ft)	POW12R1	POW12R2	PAW12R1	PAW12R2
150cm (5ft)	POW15R1	POW15R2	PAW15R1	PAW15R2
180cm (6ft)	POW18R1	POW18R2	PAW18R1	PAW18R2
210cm (7ft)	POW21R1	POW21R2	PAW21R1	PAW21R2
240cm (8ft)	POW24R1	POW24R2	PAW24R1	PAW24R2

Handrail Colour Options

Add suffix to above codes for colour option if required

Colours		Code Suffix
Original Graphite		No suffix
Original Graphite with Safety Yellow handrails		-Y
Garden Green		-G
Heritage Brown		-B

Permaramp-Quickrails dimensions:

- Add 12.5cm (5ins) to the ramp width each side the Quickrails are added
- The Quickrail length is the same as the ramp length and can be adjusted if required
- Handrail height is 92cm (36ins)
- The handrails uprights are 45cm (18ins) from each end of the ramp

Ask for more information about: Retrofitting **permaramp®-quickrails™**
Options for increased capacity • 'Portable' type applications for events, etc.

Accessories

permaramp®-quickrails™

No-Tool Kit

Code (for one side): **PRNTK**



permaramp®-entry flap

Fits standard width **permaramp®-adjust**

Code: **PE5** - See page 17



Welcome®-Giant Half-Step

Until now, the traditional half-step simply offered an ambulant user shallower step access.

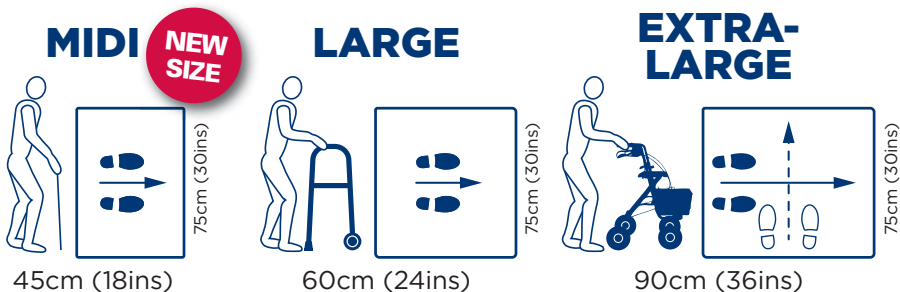
The proven and now **ALL-NEW** Welcome®-Giant Half-Step reaches a new level!

“ Making steps safe and manageable for ambulant, rollator and walking-frame users! ”

**ALL
NEW**




3 sizes...





3 challenges solved...

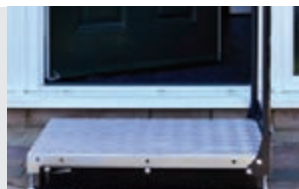
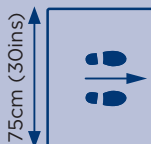
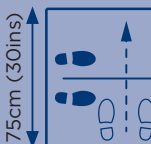
1. Safe and Manageable Step Heights








2. Eliminate the Trip Risk of Raised Threshold Obstacles

3. Instant, Safe, Stand-Alone Handrail Solution



Specifications	MIDI	LARGE	EXTRA-LARGE
 <p>Giant Half-Step With 75cm anti-trip toe-plate and rubber and fix-down feet; Height range: 8 - 15cm (3.1 - 6ins)</p>	 <p>45cm (18ins)</p>	 <p>60cm (24ins)</p>	 <p>90cm (36ins)</p>
Ambulant	✓	✓	✓
Limited Knee/Leg Mobility	-	✓	✓
Walking-Frame	-	✓	✓
Rollator	-	-	✓
Limited Space	✓	✓	-
Outward-Opening Doors	-	-	✓
Capacity	400kg (63st)	400kg (63st)	400kg (63st)

***Important note:** If using height-extension kit, 400kg capacity requires all feet to be securely fixed to the ground and requires connected handrails both sides. Capacity is otherwise 300kg (47st).

Accessories	MIDI	LARGE	EXTRA-LARGE
 <p>Handrail With integral safety kerbs. Fit to any side or opposite sides. Universal height: 90cm (36ins)</p>	WHS45 45cm (18ins)	WHS60 60cm (24ins)	WHS75 75cm (30ins) WHS90 90cm (36ins)
 <p>90cm Anti-Trip Toe-Plate Provides additional trip protection on the long (90cm) Extra-Large Half-Step sides</p>			WHSTP90
 <p>Doorplate Closes the gap between the Giant Half-Step and the door frame; self-supporting for uPVC frames</p>		WHSDP75 75cm (30ins)	
 <p>Height-Extension Kit Extend any Giant Half-Step up to 20cm (8ins) high. Useful for 'building over' existing steps. Important note: If using height-extension kit, 400kg capacity requires all feet to be securely fixed to the ground and requires connected handrails both sides. Capacity is otherwise 300kg (47st).</p>		WHSLXK	
 <p>Handrail Joiner Kit Join handrails of any two or more adjacent half-steps (2 pieces)</p>		WHSRJK	

KITS FOR EASY ORDERING	MIDI	LARGE	EXTRA-LARGE
 <p>Giant Half-Step Kit Giant Half-Step Kit with Handrail. Height range: 8 - 15cm (3.1 - 6ins)</p>	WHSMR Includes: WHSM WHS45	WHS LR Includes: WHS L WHS60	WHSX LR Includes: WHSX L WHS90
 <p>Double Giant Half-Step Kit Double Giant Half-Step kit with handrails; height up to 20cm (8ins) (includes height-extension and handrail joining kits)</p>	WHSM2R Includes: 2 x WHSM 2 x WHS45 WHSRJK WHS LKX	WHS L2R Includes: 2 x WHS L 2 x WHS60 WHSRJK WHS LKX	WHSX L2R Includes: 2 x WHSX L 2 x WHS90 WHSRJK WHS LKX

All Welcome®-Giant Half-Steps are supplied with our easy-to-follow instruction manual. For further application and assembly support, call us on 020 8275 0375

rampCENTRE® WELCOME®



Residential modular ramp & step systems

This residential ramp system has been carefully developed to suit your budget, needs and timescales. It is incredibly good value, whilst being very versatile.

Attractive • Excellent value • Easy to install • Durable • Totally modular

- Perfect for independent access
- Used when permanent access is required
- Handrail options available
- Easy to remove or modify when needs change, avoiding disruptive demolition work

Attractive modular residential ramp solution that suits your budget, needs and timescales

Your benefits:

- CAD drawings available for any variation on request
- Complies fully with Building Regulations Document M for dwellings



STRENGTH

Built to last, the system has a huge 400kg capacity.



HANDRAILS

Powder-coated pre-galvanised steel provides a durable yet 'warm to the touch' handrail.



ADJUSTABLE LEGS

Height-adjustable legs can be adjusted to suit ground variations and a huge range of step heights.



PLATFORM

2 platform sizes available. If a greater size is required, simply join them together.



DOORPLATE

Standard 75cm wide doorplate supplied with every system – self-supporting and wide versions also available (see page 50)



Fast one-person install



Small-van friendly



100% modular



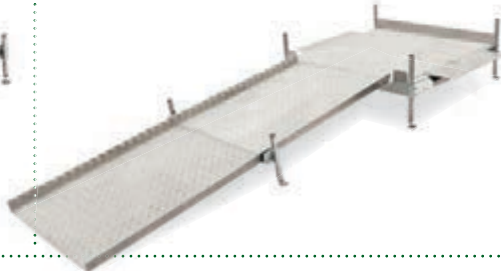
Kit A – Ramp Only

Often used to provide access to an existing platform or step.



Kit B – Ramp and Platform

A very practical installation layout, complying with the requirements for a resting platform at the top of a ramp. The 120cm square platform can be used in both straight-through and turning configurations.



Handrail Options



**No
Handrail**



**Single Handrail One
Side**



**Single Handrail
Both Sides**



**Double Handrail
One Side**



**Double Handrail
Both Sides**

WELCOME® - colours

The Welcome ramp systems are now available in a carefully selected range of colours – from stock, and at no extra cost!

**Original Graphite with
Safety Yellow handrails**



Garden Green



Heritage Brown



Colour Options



Original Graphite

The standard Welcome colour is a carefully selected and universally appealing soft dark grey.



Safety Yellow

A carefully selected rich yellow, perfect for visual safety needs. (handrails only).



Garden Green

A traditional deep green, perfect for garden environments.



Heritage Brown

A dark, rich brown to perfectly complement traditional and heritage buildings.

Choose colours to match in or contrast with every building as required!



WELCOME® - patio platform™

Suited to a wide variety of scenarios including patio, French and extra-wide doors



WELCOME® - accessories

Self-supporting doorplate for 122cm platform sides

No more cracked uPVC door thresholds!



Steps and Step Handrails

Easy access for all users.



Soft-ground leg bearer

Stop legs from sinking in soft ground.



Bow-end

A great finishing touch with safety benefits.



Welcome Platform Gate

Available in Original Graphite colour only.



Welcome Platform Door Entry Flap

COMING
SOON

Code:
WPDEF

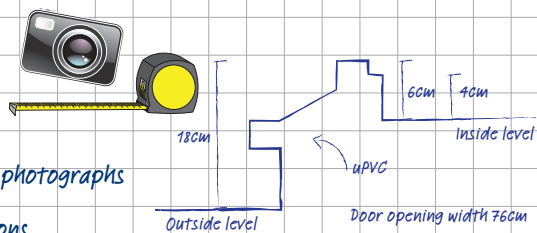


Installed Ramp Gallery



TOP TIPS!

1. Refer to Ramp Specifier Guide (p2/3)
2. Always take close-up and stand-back photographs
3. Make a rough sketch and add dimensions



rampCENTRE®

by

Enable Access

- Established 1998
- Multiple patented products
- Largest range available from stock
- Exported worldwide



“Leading Ramp and Access innovation”



Helping you
“Get Ramps Right!”

For more details contact:

Enable Access, AccessCentre, Aerolight, ExitMaster, RampCentre and Metro Mobility are registered trademarks of Enable Access. Welcome, Rollout-Trackway, Doorline, Ultralight, Permaramp, and Swift Ramp System are trademarks of Enable Access. These and all other product names, codes and images must only be used in reference to products supplied by Enable Access. We undertake a policy of continual review and research, which may lead to changes in product appearance and specification. Information is given for guidance only; customers take ultimate responsibility for purchasing decisions. Aerolight ramps are designed and optimised for disability/mobility type use - for other applications, please speak to a RampCentre product advisor. All dimensions, weights and conversions are approximate for your easy reference only and specific risk assessments (including all manual handling assessment and product use issues) should be carried out prior to the use of the equipment supplied. Stated capacities apply for combined weight of mobility devices and/or passenger/attendant. For technical help and details of your nearest stockist phone 020 8275 0375. © Copyright Enable Access.

DE-ARAG10/21

