



**FormAlign**

**New Backrest  
Range  
for Active  
Users**



**ignite**



*maximising  
postural  
efficiency*

***Stabilising one's own posture  
Constantly consumes energy***

*When stability & balance  
is improved through accurate  
trunk support, more focus is  
redirected to the task at hand.*

*From its comfortable pressure redistributing foam inner to its easy-to-use U2 compact quick release hardware & sporty appeal the light weight Dreamline Ignite Backrest System is a pleasure to use.*



*user  
focused*




# 4 pillars of design

## 1 / USER COMFORT

Utilising Dreamline's Bubble Contour design for comfortable immersion with breathability channels. A central channel is included for spinal relief. The foam pad is overlaid with a pressure redistributing viscoelastic "Sunmate" foam for extra long-term comfort.

FLO2™ engineered 4-way ultra-stretch fabrics with moisture wicking properties enable the user to fully utilise the comfort of the pad.


The pad & cover lock to the shell securely for trouble free riding. The cover is removable & machine washable.

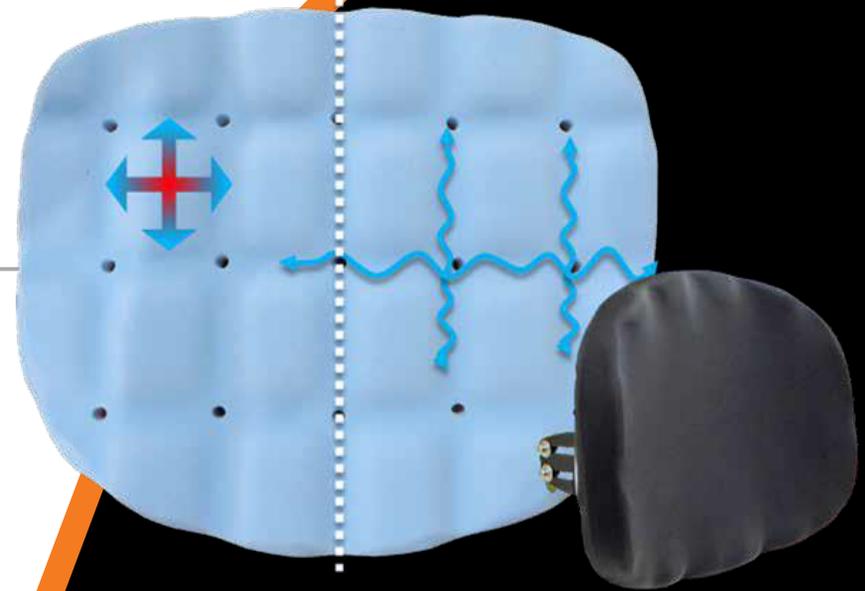

## 2 / CONTOUR PROFILE

Specifically designed with an open curvature shape with a more closed lateral angle to encourage mid-line sitting posture.

Slim 1" wide laterals maximise free arm movement & ensure a compact user-to-wheelchair fit.

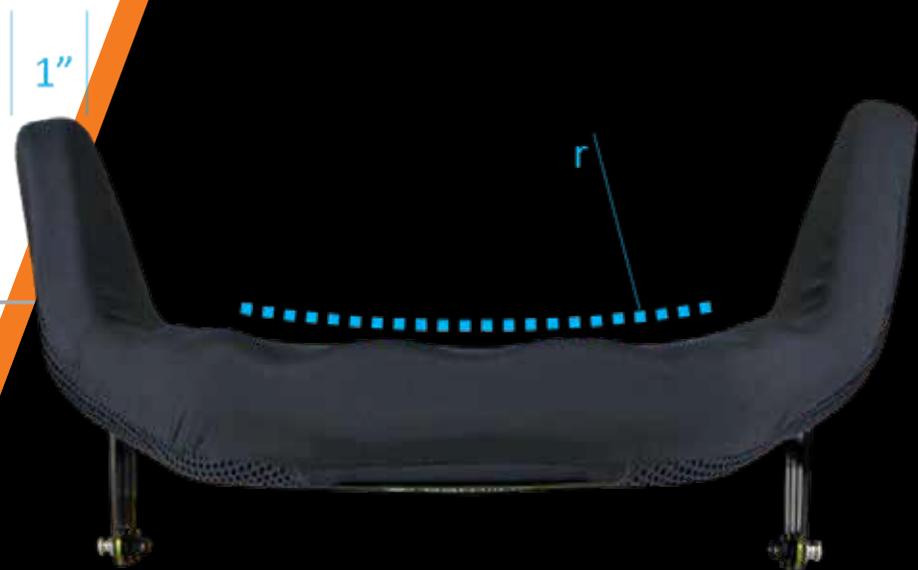



### 3 / POSITIONING

Ignite backrests have a series of positioning pads supplied in the standard kit including a Lumbar pad & a pair of wide wedges & narrow wedges. The added ability to include Dreamline's BAC Backrest Advanced Contouring system to any Ignite backrest ensures users with asymmetrical postures will be well supported. Truly reversible with an offset lateral design ensuring the lateral height best suits the user.



### 4 / ADJUSTABILITY

The proven U2 quick-release hardware is effortless & quick to operate. The U2 Compact clamps attach to tubing diameters from 19mm (3/4") to 25mm (1"). Enormous width adjustability (4") ensures the right fitting backrest is selected for the user & ensures it will fit the wheelchair frame. Up to 60deg. of angle adjustment (+/- 30deg.) & 5" of depth adjustment takes the backrests ability beyond the need.

# *easy selection*

2 types

CARBON



ALUMINIUM



3 contour  
depths






## 3 heights

1 / SHOULDER

2 / THORACIC


3 / LUMBAR

# Ignite

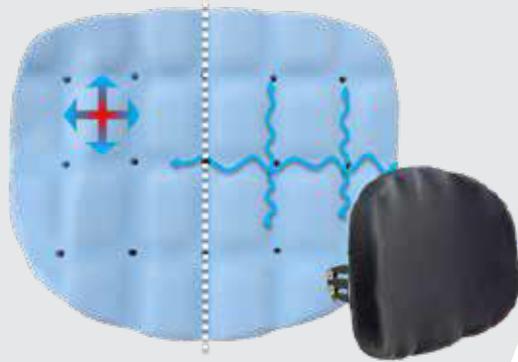


## FLO<sub>2</sub><sup>TM</sup> ENGINEERED FABRIC

FLO<sub>2</sub><sup>TM</sup> is an advanced high-performance fabric designed for medical and high energy sporting users

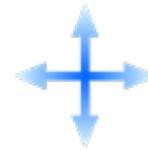

Extreme  
4-Way Stretch

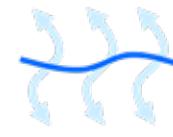

Air  
Exchange



Moisture  
Wicking



Durable



Fire  
Retardant



*tested for safety*

### **The Dreamline Ignite Backrests U2**

*compact hardware is fully crash tested to meet the standards of WC20 & ISO7176-19 for forward facing transportation with an occupant in a motor vehicle.*

*When occupied in a motor vehicle, Dreamline backrests must only be used on wheelchair frames which provide four anchorage points & pelvic belt anchor points which comply with the requirements of RESNA WC-4:2012, Section 19.*

# Sizing chart

Easy to use Prescription Form available on request

## ACTIVE + MID CONTOUR SIZE OPTIONS

		WIDTH						
		14"	15"	16"	17"	18"	19"	20"
HEIGHT	7"							
	8"							
	9"							
	10"							
	11"							
	12"							
	13"							
	14"							
	15"							
	16"							
	17"							
	18"							
	19"							
	20"							

## DEEP CONTOUR SIZE OPTIONS

		WIDTH						
		14"	15"	16"	17"	18"	19"	20"
HEIGHT	7"							
	8"							
	9"							
	10"							
	11"							
	12"							
	13"							
	14"							
	15"							
	16"							
	17"							
	18"							
	19"							
	20"							

*active / mid / deep contours*

# *empowering performance*

*with the right backrest,  
postural efficiency  
& sitting ability can  
increase, helping to  
improve overall  
performance*



Ignite



 01473 467777

 [sales@formalgroup.co.uk](mailto:sales@formalgroup.co.uk)

 Unit 7 Liberator House Clopton  
Park Woodbridge IP13 6QT

*Making Your Day!*