

# Klaxon Mobility



the art of moving



# Klick Electric Standard



Model	Klick Electric Standard
Material	Aluminium
Motor	36V x 370W brushless
Wheel	14" smooth tread tyre
Display	Multifunction LCD, 5 speed settings, walking pace
Battery	Lithium 36V x 5.8 Ah (weight 1.3 kg)
Handlebar	Aluminium, adjustable in height, depth and angle
Brakes	Double Brake system 160 mm with Electronic Brake System
Size	W 50 cm x H 90 cm x D 40 cm
Connection	Klick Linking System Klaxon® (worldwide patented)
Weight	8.0 kg (battery excluded)
Range	About 25 km (speed setting 1, level ground)
Included accessories	Battery 5.8 Ah with battery charger, reverse gear, LCD display, bell, foldable kickstand, splashguard



Linking System Klaxon® - Worldwide Patent

## EBS - Electronic Brake System

One touch of the brake lever is all it takes to stop the Klick without having to use the mechanical brake. During the EBS braking the energy produced will be recovered and used to recharge the battery. The mechanical brake can still be activated by applying a greater force on the brake lever, the two braking system will then merge together in a single braking action.

## "Walking Pace" - Speed Setting

Ideal for your relaxing strolls, the "walking pace" function will make your Klick maintain a constant speed of 4,5 km/h.





# Klick Electric Tetra



Model	Klick Electric Tetra
Material	Aluminium
Motor	36V x 370W brushless
Wheel	14" smooth tread tyre
Display	Multifunction LCD, 5 speed settings, walking pace
Battery	Lithium 36V x 5.8 Ah (weight 1.3 kg)
Handlebar	Aluminium, adjustable in height, depth and angle, Comfort Drive System (accelerator and brake activated through push and pull movement)
Brakes	Double Brake system 160 mm
Size	W 50 cm x H 95 cm x D 40 cm
Connection	Klick Linking System Klaxon® (worldwide patented)
Weight	8.6 kg (battery excluded)
Range	About 25 km (speed setting 1, level ground)
Included accessories	Comfort Drive System, battery 5.8 Ah with battery charger, reverse gear, LCD display, bell, foldable kickstand, splashguard



Linking System Klaxon® - Worldwide Patent

## "Walking Pace" - Speed Setting

Ideal for your relaxing strolls, the "walking pace" function will make your Klick maintain a constant speed of 4,5 km/h.





# Klick Race Standard



Model	Klick Race Standard
Material	Aluminium
Motor	48V x 1000W brushless, all-wheel
Wheel	14" smooth tread tyre
Display	Multifunction LCD, 5 speed settings, cruise control, walking pace
Battery	Lithium 48V X 11 Ah (weight 3.0 kg)
Handlebar	Aluminium, adjustable in height, depth and angle
Brakes	Double Brake system 160 mm with Electronic Brake System
Size	W 50 cm x H 90 cm x D 40 cm
Connection	Klick Linking System Klaxon® (worldwide patented)
Weight	11.4 kg (battery excluded)
Range	About 50 km (speed setting 1, ground level)
Included accessories	Battery 11 Ah with battery charger, reverse gear, LCD display, bell, foldable kickstand, splashguard



Linking System Klaxon® - Worldwide Patent

## EBS - Electronic Brake System

One touch of the brake lever is all it takes to stop the Klick without having to use the mechanical brake. During the EBS braking the energy produced will be recovered and used to recharge the battery. The mechanical brake can still be activated by applying a greater force on the brake lever, the two braking system will then merge together in a single braking action.

## "Walking Pace" - Speed Setting

Ideal for your relaxing strolls, the "walking pace" function will make your Klick maintain a constant speed of 4,5 km/h.






# Klick Race Tetra



Model	Klick Race Tetra
Material	Aluminium
Motor	48V x 1000W brushless, all-wheel
Wheel	14" smooth tread tyre
Display	Multifunction LCD, 5 speed settings, walking pace
Battery	Lithium 48V X 11 Ah (weight 3.0 kg)
Handlebar	Aluminium, adjustable in height, depth and angle, Comfort Drive System (accelerator and brake activated through push and pull movement)
Brakes	Double Brake system 160 mm
Size	W 50 cm x H 95 cm x D 40 cm
Connection	Klick Linking System Klaxon® (worldwide patented)
Weight	12.0 kg (battery excluded)
Range	About 50 km (speed setting 1, ground level)
Included accessories	Comfort Drive System, battery 11 Ah with battery charger, reverse gear, LCD display, bell, foldable kickstand, splashguard

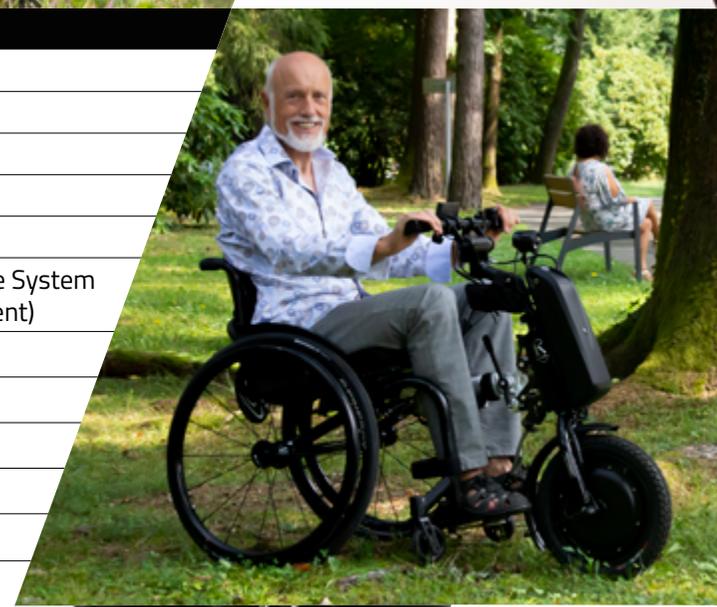

Linking System Klaxon® - Worldwide Patent

## "Walking Pace" - Speed Setting

Ideal for your relaxing strolls, the "walking pace" function will make your Klick maintain a constant speed of 4,5 km/h.





# Klick Mini Standard



Model	Klick Mini Standard
Material	Aluminium (with standard carbon cover)
Motor	36V x 370W brushless, all-wheel
Wheel	8" smooth tread tyre
Display	Multifunction LCD, 5 speed settings, walking pace
Battery	Lithium 36V x 5.8 Ah (weight 1.3 kg)
Handlebar	Foldable handlebar and steer
Brakes	Single brake 140 mm
Maximum Dimension	(when open) W 45cm x H 90cm x D 25cm
Minimum Dimension	(when closed) W 25cm x H 65cm x D 25cm
Connection	Klick Linking System Klaxon® (worldwide patented)
Weight	7.7 kg (battery excluded)
Range	About 25 km (speed setting 1, level ground)
Included accessories	Battery 5.8 Ah with battery charger, LCD display, foldable kickstand



Linking System Klaxon® - Worldwide Patent

## "Walking Pace" - Speed Setting

Ideal for your relaxing strolls, the "walking pace" function will make your Klick maintain a constant speed of 4,5 km/h.






klaxon-klick.com  
KLAXON

# Klick Carbon Standard



Model	Klick Carbon Standard
Material	Aluminium (with Carbon Cover, Led light integrated)
Motor	36V x 370W brushless, all-wheel
Wheel	8" smooth tread tyre
Display	Multifunction LCD, 5 speed settings, walking pace
Battery	Lithium 36V x 5.8 Ah (weight 1.3 kg)
Handlebar	Foldable handlebar and steer
Brakes	Single brake 140 mm
Maximum Dimension	(when open) W 45cm x H 90cm x D 25cm
Minimum Dimension	(when closed) W 25cm x H 65cm x D 25cm
Connection	Klick Linking System Klaxon® (worldwide patented)
Weight	7.9 kg (battery excluded)
Range	About 25 km (speed setting 1, level ground)
Included accessories	Battery 5.8 Ah with battery charger, LCD display, Led light integrated, foldable kickstand

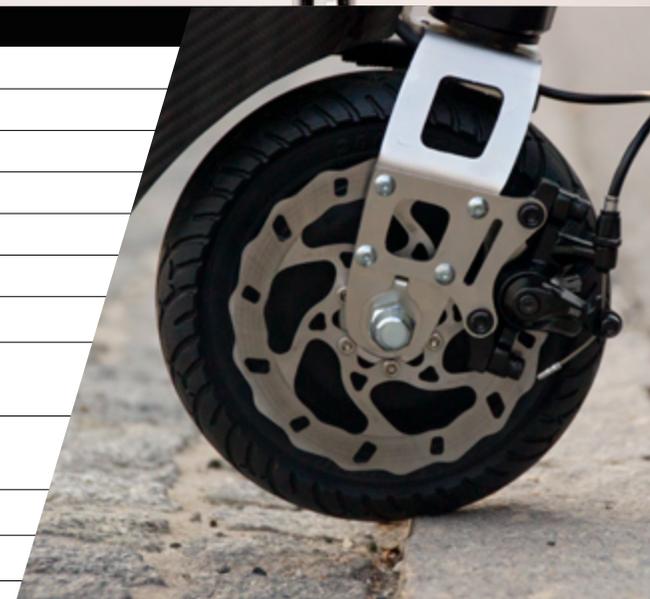

Linking System Klaxon® - Worldwide Patent

## "Walking Pace" - Speed Setting

Ideal for your relaxing strolls, the "walking pace" function will make your Klick maintain a constant speed of 4,5 km/h.





# Klick Monster Standard



Model	Klick Monster Standard
Material	Aluminium
Motor	48V x 1000W brushless
Wheel	20" tyre XL
Display	Multifunction LCD, 5 speed settings, cruise control, walking pace
Battery	Lithium 48V x 11 Ah (weight 3.0 kg)
Handlebar	Aluminium, adjustable in height, depth and angle
Brakes	Double Brake system 160 mm with Electronic Brake System
Size	W 50 cm x H 90 cm x D 50 cm
Connection	Klick Linking System Klaxon® (worldwide patented)
Weight	15.0 kg (battery excluded)
Range	About 40 km (speed setting 1, ground level)
Included accessories	Battery 11 Ah with battery charger, reverse gear, LCD display, bell, foldable kickstand, splashguard



Linking System Klaxon® - Worldwide Patent

## EBS - Electronic Brake System

One touch of the brake lever is all it takes to stop the Klick without having to use the mechanical brake. During the EBS braking the energy produced will be recovered and used to recharge the battery. The mechanical brake can still be activated by applying a greater force on the brake lever, the two braking system will then merge together in a single braking action.

## "Walking Pace" - Speed Setting

Ideal for your relaxing strolls, the "walking pace" function will make your Klick maintain a constant speed of 4,5 km/h.





# Klick Hybrid Standard

Model	Klick Hybrid Standard
Material	Aluminium
Motor	48V x 500W brushless
Wheel	20" smooth tread tyre
Display	Multifunction LCD, 5 push-assistance selection
Battery	Lithium 48V x 11 Ah (weight 3.0 kg)
Handlebar	Aluminium adjustable in height
Brakes	Double Brake system 160 mm
Size	W 60 cm x H 125 cm x D 50 cm
Connection	Klick Linking System Klaxon® (worldwide patented)
Weight	18.0 kg (battery excluded)
Range	About 50 km (in electric mode, speed setting 1, level ground)
Gear	9 x 1
Included accessories	Battery 11 Ah with battery charger, LCD display, USB-port phone charger, foldable kickstand, splashguard



Linking System Klaxon® - Worldwide Patent

### Quick and easy to attach

The innovative linking system allows for your favourite KLICK handbike to be connected and disconnected within seconds without any personal assistance, so you can explore new horizons and move independently.



### Universal connection hub

The KLICK series is suitable for all manual wheelchair types, including foldable wheelchairs, thanks to the universal connection hub. Experience a new level of freedom without compromising the manoeuvrability and comfort of your manual wheelchair.



### Suitable for different settings

Your KLICK can take on any environments, providing you a great terrain coverage, whether you want to climb steep roads, go on an off-road adventure or manoeuvre easily through lively shopping streets or small indoor spaces.



*"I want to make memories...*

*all over the world..."*



### Safe and secure

With our high level of quality and performance we want to ensure that your KLICK takes you where you want to go, so you can move independently and engage in a more active lifestyle.



### Light and easy to handle

Experience a new level of freedom, as your KLICK will help you overcome everyday obstacles, climb steep roads and manage difficult terrains with ease. It is easily manoeuvrable, allowing you to turn in small places.



### Compact and suitable for travel

Explore new horizons, go on an adventurous trip or just cruise the streets for the pure joy of it! You will enjoy it!



1



19



18



16



15



13



17



14



12

## Main Accessories and Optionals



2

1 Mini Flybag - FLYBAG01

2 Standard Flybag - FLYBAG00

3 Race Battery (48V x 11Ah) - CM-1-11AH  
Fly-Version (280 Wh) - CM-1-5AH

4 Slim Battery (36V x 5.8Ah) - CC-1-5AH  
Fly-Version (90 Wh) - CC-1-2,5AH SLIM

5 Power Battery (36V x 11Ah) - CC-1-11AH  
Fly-Version (90 Wh) - CC-1-2,5AH PLUS

6 Mini Battery (36V x 5.8Ah) - MK-1-5AH  
Fly-Version (90 Wh) - MK-1-2,5AH

7 Hybrid Battery (48V x 11Ah) - CH-1-11AH

8 Quick Release Adjustable Handlebar Stem  
CC-23-22-24-0001

9 Klaxon Gloves - GLOVES

10 Addictional Traction Weight - CC-93-0000

11 Phone Holder - PHONE-HOLDER

12 Standard Light - LED-WHITE-0000

13 Front Light Monster/Race - CC-83-KIT

14 Handlebar Extension Std - CC-25-0005  
Handlebar Extension Tetra - CC-25-0004

15 Race Battery Charger - CC-88-0001

16 Power Battery Charger - CC-88-0000  
Hybrid Battery Charger - K500050

17 Slim Battery Charger - CC-89-0000  
Mini Battery Charger - MK-88-0000

18 X-Large Wheels Set  
FATWHEELS

19 14" Tyres - CC-3-1400



3



4



5



6



7



8



10



9

KLAXON

# Linking System Klaxon®

Worldwide Patented



The innovative linking system connects to any hand bike of the KLICK product range in seconds, allowing you to switch between your favourite KLICK hand bikes with ease. The universal connection hub adapts to all manual wheelchair types, including foldable wheelchairs. Experience a new level of freedom without compromising the manoeuvrability and comfort of your manual wheelchair.

**LINE-UP. ENGAGE. KLICK.**

QUICK  
and easy to attach



LINKS  
in seconds



UNIVERSAL  
connection hub



[www.klaxon-klick.com](http://www.klaxon-klick.com)

[info@klaxon-klick.com](mailto:info@klaxon-klick.com)



Facebook:  
[@KlaxonKlickWorld](https://www.facebook.com/KlaxonKlickWorld)

Instagram:  
[@KlaxonKlick](https://www.instagram.com/KlaxonKlick)