

OPTIMUS 2 2.322


USER PROFILE

- ✓ For people using a manual wheelchair or walking aid indoors and who require more support outdoors

HIGH-PERFORMANCE OUTDOOR POWER WHEELCHAIR (WITH FRONT-WHEEL DRIVE)

More powerful batteries and a top speed of up to 15 km/h make the all-terrain power wheelchair an experience.

- Ideal for long distances out of doors: travel up to 111 km on a single charge with the 6 km/h version.*
- Optimised chassis geometry enables best-possible ride quality
- Single-wheel suspension and high climbing ability are ideal for high kerbs, higher slopes and challenging and rough ground
- Direct steering combined with a travel speed of 6, 10 or 15 km/h for high-level safety and reliability outdoors
- Larger batteries available for long range and outstanding climbing ability
- Extensive range of seats for optimally adapted fitting

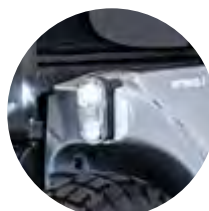

LED front lighting



Obstacle climbing ability up to 110 mm



Electrical seat and back angle adjustment also for more complex fittings (optional)

R-Net electronics as standard

Puncture-proof tyres with studded tread (optional)

Drum brakes as standard



User weight up to 150 kg



Very powerful motor 950 W



Very high max. range between 60 – 111 km (depending on battery size and speed)



* Tested according to: DIN EN 12184:2022 'Electrically powered wheelchairs, scooters and their chargers – Requirements and test methods Chapter 8.1.9 Range in connection with ISO 7176-4:2008 'Energy consumption of electrically powered wheelchairs and scooters for determining theoretical range'

DETAILS

OPTIMUS 2 2.322

MORE FREEDOM WITH THE ENERGY- SAVING RANGE WONDER

Thanks to optimisations to the OPTIMUS 2, the 15 km/h version is now even more powerful with a range of up to 77 km. For better visibility on dull days, bright LED lighting has also been installed. Quality has also been significantly improved with the new robust iCHAIR MEYLIFE legrests. High climbing ability of 18% and low power consumption during the charging process are further important advantages.

New 15 km/h software

Improved drive stability and range (up to 77 km in the 15 km/h version with 125 AH (20h) battery),

**approved for
18% gradient**



LED lighting

Bright front and rear lighting and new white front safety reflector



Integrated legrest

Detachable, robust legrests



Mechanical
height adjustment



Electrical
height adjustment