

Q700M SERIES

Premium comfort, from sit to stand.



Q Series
Powerchairs



Q700 M Pro



Q700 M Advanced Pro



Q700-UP M Advanced

Introducing the **QUICKIE Q700 M series** - enjoy superior outdoor performance and enhanced comfort without sacrificing the compact base for tight indoor manoeuvrability and intuitive driving that you love from a mid-wheel drive powerchair.

Empower your mobility experience without the need to compromise on either comfort or breaking the bank with our selection of tailored seating options – the **SEDEO PRO**, the **SEDEO ADVANCED PRO**, or **SEDEO UP**.



Scan the QR code with
your phone's camera to
book a demo.

Q700M

SERIES OVERVIEW

PERSONALISED COMFORT, YOUR WAY.

Enjoy the freedom to select the perfect fit for your individual lifestyle needs within the Q700 M series, where comfort meets adaptability.



Q700 M PRO

With essential power features like lift, tilt, anti-shear recline, and power leg rests, empowering you to effortlessly reposition at the touch of a button, ensuring comfort tailored to your needs.



Q700 M ADVANCED PRO

Designed for superior functionality. Boasting the same features as the Q700 M PRO but with added benefits like lay-flat and partial stand, promoting excellent posture even during transitions. Providing unparalleled comfort and support.



Q700-UP M ADVANCED

With all the features of the Q700 M ADVANCED PRO but loaded with cutting-edge standing technology, this model brings premium lifestyle and postural features to you at an affordable price point. Experience the next level of mobility with unmatched value.





DRIVE BASE

TRUE MID-WHEEL DRIVE TECHNOLOGY.

The Q700 M series shares one drive base that sets a new standard both inside and out. This series seamlessly blends unparalleled indoor manoeuvrability with exceptional outdoor performance, ensuring your comfort, safety, and stability on the most challenging terrains.

INDOOR PERFORMANCE

SUPERIOR INDOOR AGILITY.

Staying in? Make yourself at home with unparalleled ease of movement. The Q700 M series offers superior indoor manoeuvrability and intuitive driving, allowing you to effortlessly navigate tight, compact spaces.





MOVE ANYWHERE

Tired of limitations? Whether it's a challenging living room, a busy kitchen, or a narrow doorway, the Q700 M's compact **620 mm** base seamlessly fits into your life, ensuring you move freely without constraints.



COMPACT 1100 MM TURNING CIRCLE

With its **TRUE** mid-wheel drive technology, the Q700 M effortlessly turns on the spot, boasting a compact **1100 mm** turning circle. In contrast, RWD and FWD powerchairs need to 'swing' around on their front or rear drive wheels, resulting in a wider turning circle.

OUTDOOR PERFORMANCE

BUILT FOR THE OUTDOORS.

The ideal companion for conquering the outdoors – experience extended travel range with superior drive performance and feel confidently secure even on the most challenging landscapes and while navigating obstacles.

ENHANCED TRACTION

The Q700 M's mid-wheel drive technology strategically places the drive wheel directly beneath your seated position, ensuring more weight is on the drive wheel. This design delivers superior traction compared to front or rear-wheel drive bases, providing you enhanced stability and control.





UNLEASH MAXIMUM FREEDOM AND POWER

Experience the thrill of power and speed with the 4-pole high-torque motors, placing control at your fingertips. Standard options of 6, 8 or 10 for ultimate performance. Choose between 60 Ah or 80 Ah batteries to maximize your range, reaching up to **40 km** and unlocking unparalleled outdoor performance.

GO ANYWHERE WHEELS

The Q700 M series is equipped with **14"** drive wheels which offer unmatched performance. Engineered for off-road driving, these wheels conquer the common challenge of getting stuck where others would.





ADVANCED SUSPENSION

TAILORED SUSPENSION FOR YOUR COMFORT.

Say goodbye to discomfort during your drive. The Q700 M ensures an unprecedented level of comfort with independent articulation of each wheel. Coupled with its unique suspension system, navigating obstacles and traversing rough terrain has never been smoother, providing you with a truly comfortable and enjoyable experience.

ALWAYS IN CONTROL

QUICKIE's anti-pitch technology eliminates pitching when decelerating, meaning you always feel safe and in control. With all your weight located directly above the drive wheel to maximise traction and stability, you've never been in a safer driving position.

ADVANCED ALL-WHEEL SUSPENSION

Experience unparalleled comfort on every journey as the Q700 M's all-wheel suspension effortlessly glides over obstacles and bumps, ensuring you barely feel a thing. Optimising ride comfort, the independent movement of all six wheels empowers the Q700 M to conquer any terrain, providing a journey where comfort meets security.



CLIMB WITH COMFORT

Conquer obstacles up to **75 mm** effortlessly, and revel in a smoother, more responsive ride with revolutionary suspension. The independent movement of all six wheels, including both front and rear, gives the Q700 M the ability to crawl up and down steep obstacles while ensuring all six wheels stay grounded. Maximising comfort and traction, even on gradients.



C-ME

ELEVATE AND DRIVE


Say goodbye to waiting around or driving at reduced speed when elevated. With C-ME*, you can raise up to 300 mm and drive up to 4.8 kph all at the same time! Now you can keep up with your friends at their pace and be at eye-level all with a touch of a button (with no need to stop).



Scan the QR code with your phone to learn more about C-Me Technology with our informative video

RAISING COMFORT WHEN DRIVING ELEVATED

Unlike other powerchairs, our suspension doesn't lock-out when using C-ME. So you get the same smooth and comfortable ride over uneven terrain, even when elevated! Plus, anti-pitch technology will keep you safe and stable on inclines too.

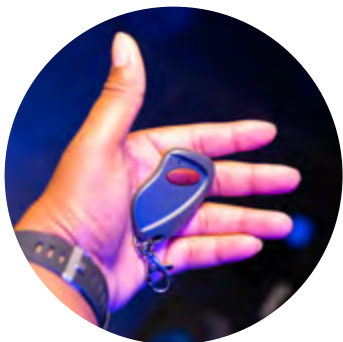

STOP SWITCH FOR ADDED SAFETY

There's safety in stopping. Our added remote stop switch allows an attendant to stop a powered wheelchair via a remote from up to **30 meters** away. Operating distance dependent on the environment.



ELEVATE WHILST DRIVING WITH ASSIGNABLE BUTTONS

Fed up of complex joystick movements or having to stop to use your seat lift? With our **PATENTED** assignable buttons, you can change any button (or assign a second, long-press function) to activate C-ME (even when driving).











DISTINGUISHED STYLE

MAKE IT YOUR OWN.

With the Q700 M series, we've taken style to the next level. With a matt-black chassis finish, hidden wiring and components combined with sharp-tailored upholstery, it's a sleek and contemporary design.

FANCY A CHANGE?

The Q700 M series also comes with a choice of 7 contemporary colour inlays to perfectly complement its striking matt-black bodywork. Simply replace the colour inlays. Styling your powered wheelchair has never been so simple.

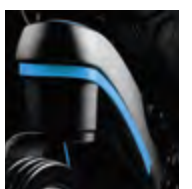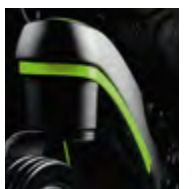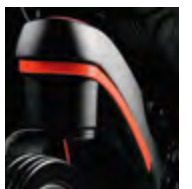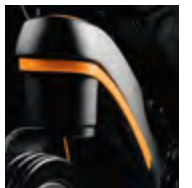

SEATING

YOUR LIFE, YOUR WAY. SEATING TAILORED TO YOU.

When it comes to selecting a powerchair, why compromise your lifestyle for clinical needs?

With three versatile seating systems, the Q700 M series ensures you don't have to! Whether you prefer a chair that tilts and lifts, one that allows you to lie flat, or the option to stand in your chair, the Q700 M series caters to your preferences.

Each seat is prescribed to your dimensions and can be adjusted as you grow. Furthermore, all options come with a comprehensive range of accessories to support your posture, maximising your daily independence.







SHARED ACROSS THE RANGE

CHOOSE YOUR COMFORT ZONE.

The Q700 M provides a choice of three seating options: SEDEO PRO, SEDEO PRO ADVANCED and SEDEO UP. Select the one that aligns seamlessly with your postural, lifestyle, and budget preferences.





TAILORED TO YOUR POSTURE AND LIFESTYLE

The Q700 M series presents an extensive array of seat accessories, empowering you to personalise your Q700 M according to your unique requirements. Choose from a variety of options, including power functions (lift and tilt, recline, legrests), armrests, armpads, headrests, and your preferred style of legrests and much more.



DESIGNED TO FIT YOU

The Q700 M series is individually prescribed to your precise measurements, guaranteeing the utmost support from your seat. Should your needs evolve, the seating system can be easily adapted, accommodating seat widths from 15" to 22" and seat depths from 16" to 22".



REFINED COMFORT SEATING

Experience the epitome of comfort with the Q700 M series. Our SEDO PRO cushion and backrest contours to the shape of your body for unparalleled immersion that must be felt to be believed. Revel in comfort with the Q700 M series.

Q700 M SEDEO PRO SEATING

SEATING ADJUSTED TO YOU.

The SEDEO PRO seating system, is meticulously crafted for a personalised experience. With a highly adjustable seating frame and a wide array of positioning and lifestyle accessories, the Q700 M SEDEO PRO can be completely tailored to provide a precise and comfortable fit just for you.



Scan the QR code with your phone to learn more about SEDEO PRO Flexible Back.




FLEXI BACKREST

Tired of uncomfortable laterals and firm backrests to support your posture? The **SEDEO FLEXI** backrest effortlessly combines adjustable full rear and side trunk support with the highest level of immersion that has the be felt to be believed.



SHAPES TO YOUR SPINE

Instead of a solid base, the flexible back cushion has 8 adjustable straps that wrap around your back. Where other backrests make limited contact with the spine, the flexible back cushion's layer of soft foam contours to you, providing complete contact, comfort and stability.

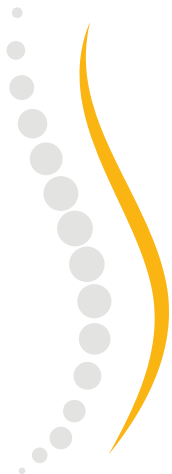

POWERED TO YOUR NEEDS

SEDEO PRO puts the power in your hands with a versatile range of functions that can be individually or collectively configured to achieve a near lie-down position:



Anti-Shear Recline - Adjust your backrest electronically through the joystick allowing you to independantly reposition yourself throughout the day. The anti-shear mechanism means there is no sliding between you and the backrest when adjustment ensuring comfort throughout use. You can also combine this with power legrests to allow you to be positioned in a near lie down position for greater comfort and pressure relief.



Lift & Tilt - Designed for relieving pressure by tilting reward. As well as offering the ability to lift the seat so you can more easily access shelves, cupboards and talk to friends eye to eye.



Powered Legrest - Adjust your legrest angle quickly through the joystick for comfort, pressure relief or oedema control when combined with tilt.

Q700 M SEDEO ADVANCED PRO SEATING

LIE-FLAT, PARTIAL STAND, AND MORE.

Elevate your seating experience with SEDEO ADVANCED PRO, featuring powered memory functions that customise seating positions to your preferences. Effortlessly access your favourite positions, including lying flat, partial standing, pressure relief, and more, ensuring comfort at your fingertips.



Scan the QR code with your phone to learn more about JAY SYNCRA.



SEATING THAT REMEMBERS YOU

Say goodbye to the frustration of manually adjusting each power function on your chair. With the ADVANCED PRO, you no longer have to. Save up to 6 tailored positions, precisely customised to your preferences - be it lying flat, partial standing, standard seating position, or pressure relief. Paired with QUICKIE's assignable buttons, you can effortlessly access these positions with just a touch of a button, making your experience easy and quick.

LIE-FLAT AND RELAX

Experience relaxation with the SEDEO ADVANCED PRO, by going in to full lay flat or relaxed position. Transitioning in and out of recline is made effortlessly comfortable with 180 mm of anti-shear, ensuring your seated posture, headrest position, and laterals are maintained throughout the repositioning process.



EMPOWER YOUR EVERYDAY TASKS WITH PARTIAL STAND

Revolutionise how you access the world with Partial Stand integrated into the ADVANCED PRO seating system. Partial Stand tilts the seat forward, raises the back, lowers the front, and positions the knees below the pelvis. This innovative function enables easy access to tables, cupboards, and other everyday tasks, promoting eye-to-eye communication and providing pressure relief for added comfort.

SEATING IN HARMONY WITH YOUR MOVEMENTS

Introducing a first-class seat meticulously designed to move with you – the patented Biometric Repositioning. Uniquely attuned to the body's movements, this seating system is among the most effective for independent repositioning throughout the day. Experience a consistently comfortable seated position, even after adjusting your seat. This is accomplished by transitioning through an intermediary tilted position, utilising gravity for a seamless realignment of your pelvis in one effortless motion.



Q700 M SEDEO PRO ADVANCED UP SEATING

RISE TO NEW HEIGHTS WITH THE Q700 UP M.

Explore life from a fresh perspective with SEDEO UP, seamlessly integrating seating system with the ability to stand at your convenience.

- 6 tailored memory positions
- Seating that mirrors the bodies movement
- Ability to change your seat position at the touch of a button
- Lie flat position, semi stand, pressure relief and more



Scan the QR code with
your phone to learn more
about JAY SYNCRA.



JAY SYNCRA

JAY SYNCRA is designed specifically to optimise your seated position whilst using powered tilt and recline. The SYNCRA seat and backrest are designed to minimise the movement of pelvis maximising your comfort and supporting a better seated position.

This is achieved through a tapered backrest, scaloped lower back shell and a lowered rear pelvic wall on the cushion and segmented fluid insert.



BENEFITS OF STANDING

With a standing angle of up to 85°, the SEDEO PRO ADVANCED UP allows you to stand whenever you want and offers many health and lifestyle benefits:

- Improved pressure relief for those who may have limited ability to independently weight shift or have to sit for long periods.
- Enhanced function and ability to reach and access your environment.
- Improved range of motion and a reduction of the risk of contractures.
- Supporting vital organ capacity, including pulmonary, bowel and bladder function and improves circulation.
- Promoting bone health and reducing the occurrence of skeletal deformities by prolonged sitting.
- Reducing abnormal muscle tone and spasticity.



SERIES OVERVIEW

WHICH WILL YOU CHOOSE?

Now that you've taken a comprehensive look at the entire Q700 M series, which one will you choose to experience firsthand?



Q700 M PRO

With essential power features like lift, tilt, anti-shear recline, and power leg rests, empowering you to effortlessly reposition at the touch of a button, ensuring comfort tailored to your needs.



Q700 M ADVANCED PRO

Designed for superior functionality. Boasting the same features as the Q700 M PRO but with added benefits like lay-flat and partial stand, promoting excellent posture even during transitions. Providing unparalleled comfort and support.



Q700-UP M ADVANCED

With all the features of the Q700 M ADVANCED PRO but loaded with cutting-edge standing technology, this model brings premium lifestyle and postural features to you at an affordable price point. Experience the next level of mobility with unmatched value.

FEATURES COMPARISON

	SEDEO PRO	SEDEO PRO ADVANCED	SEDEO PRO ADVANCED UP
Q700 M Base	✓	✓	✓
Power rear tilt	✓	✓	✓
Lift	✓	✓	✓
C-ME (3.5mph driving in full lift)	✓	✓	
SEDEO Flexi backrest	✓	✓	✓
JAY Syncra		✓	✓
Full Lie flat		✓	✓
Assignable button (change your position at the touch of a button)		✓	✓
Memory positions		✓	✓
Partial stand		✓	✓
Full stand			✓



TECHNICAL DATA

	Q700 M PRO	Q700 M ADVANCED PRO	Q700-UP M ADVANCED
Seat width	380 - 560 mm	380 - 560 mm	380 - 560 mm
Seat depth	400 - 560 mm	400 - 560 mm	400 - 540 mm
Back height	530 - 720 mm	560 - 760 mm	400 - 540 mm
Total width	622 - 660 mm	622 - 660 mm	622 - 660 mm
Total length	max. 1100 - 1190 mm	max. 1100 - 1190 mm	max. 1100 - 1190 mm
Turning circle	560 - 630 mm	560 - 630 mm	560 - 630 mm
Powered recline	90° / 150°	85° / 175°	85° / 175°
Maximum safe slope	6°	6°	6°
Maximum kerb climbing	100 mm	100 mm	80 mm
Maximum range (ISO 7176-4)			40 km (with 80 Ah batteries)
Maximum user weight	160 kg	160 kg	120 kg
Drive wheels	14"	14"	14"
Castor wheels	6"	6"	6"
Maximum speed			6, 8 and 10 kph
Battery size			60 Ah / 80 Ah (optional)
Electronics			R-Net Controller
Crash tested	Yes - approved for forward facing transport in a vehicle (ISO 7176-19). Attachment: 4 chair fixing points with 2 double straps rear & 2 single straps front. Overall rating of belt restraint accommodation according ISO 7176-19:2008 Annex D=A. No anchored belt restraints allowed. The ease of access and manoeuvrability in vehicles depends on the wheelchair size.		



IMPROVING
PEOPLE'S
LIVES®

Thorns Road, Brierley Hill, West Midlands, DY5 2LD
t: +44 (0) 845 605 6688 | e: marketing@sunmed.co.uk

sunrisemedical.co.uk

Sunrise Medical cares about your safety. Please read your owner's manual for safe operating instructions and maintenance guidelines. In line with our policy of continual product improvement, Sunrise Medical reserves the right to make changes at any time, and without notice. For the latest information, please visit www.sunrisemedical.co.uk. If you are visually impaired, this document can be viewed in PDF format at www.sunrisemedical.co.uk

