



## Bariatric Shower Commode Tilt

Bariatric Shower Commode Electric Tilt has a 32 degree tilt position and a weight capacity of 325 kg Safe Working Load. This shower commode is designed to provide a comfortable alternative showering system in instances where the prone position might cause discomfort or distress to the user. It may also be used in circumstances of limited space providing a bathing solution where other units may not have room for operation.

### **Making tasks easier and safer**

Bariatric Shower Commode Tilt provides better access for the caregiver in terms of hygiene, skin and tissue care as well as nursing or dressing feet in standing position, allowing the caregiver to work in a safe and ergonomically sound way.

### **Improving care, efficiency and communication**

This shower commode raises the efficiency of care because it enables a single caregiver to perform the whole showering procedure, as well as a wide range of everyday hygiene tasks. The high level of support provided by the shower commode results in streamlined routines that offer considerable time savings compared with traditional methods. It also enables the caregiver and the user to maintain eye contact for better communication.

### **Suitable for**

Shower Commode Tilt is suitable for patients with reduced muscle tone, pressure ulcers, spasm, involuntary movement, amputations.

- Electric tilt up to 32°
- High backrest and sturdy armrests provide optimum user support and security
- Reclined position gives full access for the caregiver to most parts of the users body.

**MAX 325 kg**

### **Product video**

Scan the code and see the product in use



With standard footrests

## Bariatric Shower Commode Tilt Specifications

Bariatric Shower Commode Tilt	
Maximum weight capacity is	325 kg
Seat depth adjustable	42 / 48 / 54 cm
Total length (lowest position)	140 cm
Total width	74 cm
Seat width	66 cm
Seat height	51-69 cm
Backrest	24-53 cm
Width of footrest	58-110 cm
Height of footrests	40-48 cm
Electric tilt	0-32°
Width between standard armrests	61 cm
Width between butterfly armrests	76 cm
Width toilet opening	Diameter 24 cm, 35 cm long
Distance from floor to seat underside	44,5 cm
Pushbar from floor	100 cm
Seat to footrest	39-42-45-48 cm
Height seat to armrest	22 cm
Weight	39 kg
Swing-away armrests	Yes
Swing-away footrests	Yes
Removable footrests	Yes
Height adjustable footrests	Yes
Headrest	Yes
4 wheels with individual locks	Yes
Control	Via the control panel
Capacity up/down fully charged battery	Approx. 20 times
Turning radius	158 cm
Frame	Powder lacquered chromium molybdenum steel
Wheels diameter	150 mm
Seat / back / armrests	Polyurethane
IP - classification	IP44
ISO - classification	ISO 09 12 03 - 23

### Ideal use for:


- Bariatric users
- Users with reduced muscle tone, pressure sores, spasm or amputations.



Headrest



Seat with lid

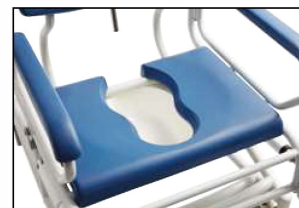

Split toilet seat

Model	Seat width cm	Seat depth cm	Order No
XXL-Rehab Shower Commode Tilt incl. standard footrests, charger, headrest and commode bucket	66	42 / 48 / 54	0155-066-000
<b>Accessories</b>			
Butterfly armrest left (add 14 cm to seat width)			0150-999-012L
Butterfly armrest right (add 14 cm to seat width)			0150-999-012R
Headrest adjustable			0156-999-020
Full size seat for shower commode with lid for bucket			0150-999-023
Commode bucket			0140-999-001
Standard footrest right			0156-999-090H
Standard footrest left			0156-999-090V
Calf support / amputation support right			0156-999-030
Calf support / amputation support left			0156-999-031



The handset gives the caregiver control over battery-powered functions