



ZIPPIE voyage

THE PERFECT BLEND
OF STYLE AND
FUNCTIONALITY



MODERN. EASY TO USE. SUPER STYLISH

The Zippie Voyage early intervention stroller combines Baby Jogger™ stroller expertise with Zippie's incredibly versatile seating system. The result is an easy-to-use stroller that makes no compromises when it comes to what's best for you and your baby.



baby
jogger™

// MODERN STYLING

The Zippie Voyage uses the popular Baby Jogger™ City Select base with its sleek styling and lightweight design. The base has been specially adapted to incorporate the Zippie Voyage seat system.

// TELESCOPIC STROLLER HANDLE

Responsive, one-hand steering that adjusts to the caregiver's height.

// TWO STEP FOLD

The frame folds into a very compact package in two fast and simple steps. Once the seat is removed then just pull upwards on both handles.



// SOFT BREATHABLE PADDING

Soft polyurethane foam enhances comfort whilst the spacer fabric is both durable and breathable to keep your baby cool and comfortable.

// LARGE STORAGE BASKET


Basket conveniently stores up to 7kg of your essential items. Low down and out of sight.

// SWIVEL FRONT WHEELS WITH SUSPENSION AND LOCK-OUT

Glide over and around obstacles while keeping baby calm, secure, and comfortable. Lock-out the front wheels for smooth tracking during long-distance strolls or easier access in and out of adapted vehicles.

ZIPPIE voyage



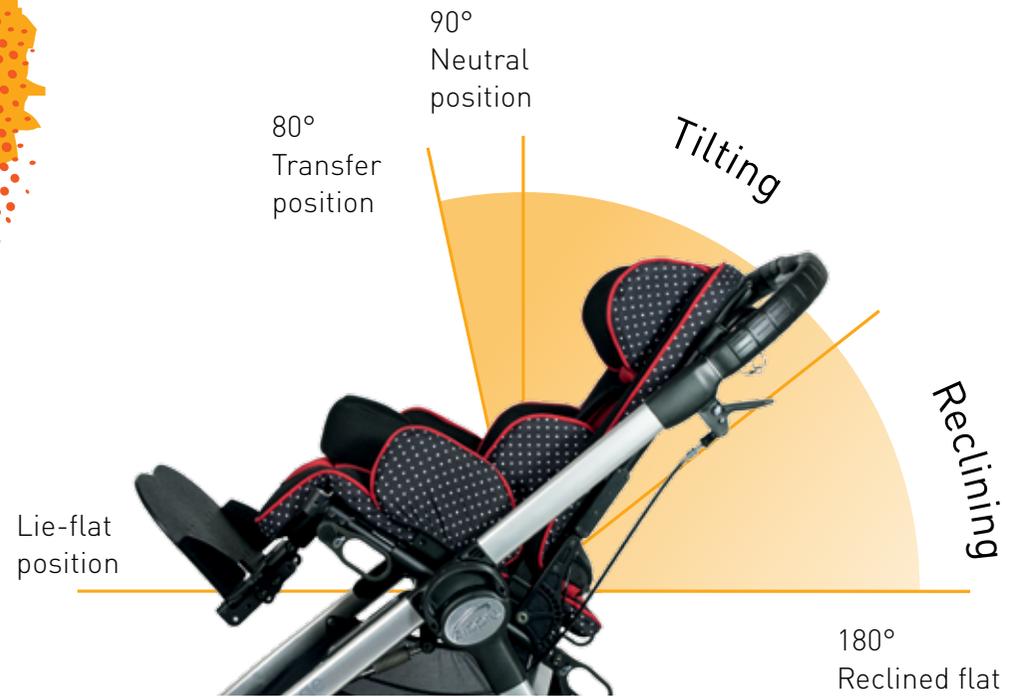

// SIMPLE TILT AND RECLINE ADJUSTMENT

Tilt and recline improves positioning and helps with feeding, digestion, respiratory function and visual orientation. Both tilt (up to 35°) and recline (up to 180°) are activated by using easy-to-reach trigger handles behind the backrest.

SUPPORTIVE SEATING FOR YOUR CHILD WITH NO COMPROMISES!

When it comes to seating, the Zippie Voyage's number one priority is positioning your child for healthy growth and development from babyhood into their school years.

// JAY™ PELVIC SUPPORT BELTS
A choice of two or four point padded positioning belts for greater pelvic control.

// POMMEL
Depth and height adjustable to promote good thigh position.

// FOOTREST
Height and angle adjustable footrest from 15.5 - 31cm, with optional foot straps and shoe holders for extra support.

// HEAD SUPPORT
Adjustable in width, height and angle to promote a good head position.

// JAY™ CHEST SUPPORT HARNESSSES
A choice of standard or dynamic chest harnesses to provide increased trunk control.

// BACKREST
Easily adjustable in height to accommodate varying trunk heights from 41 - 61cm.

// LATERAL TRUNK SUPPORT
Easily width, height and angle adjustable to ensure correct support to the trunk.

// HIP GUIDES
Width, depth and angle adjustable (from 15.5cm - 30cm) for greater pelvic and thigh support.

// SEAT PLATE
Easily adjust seat depth from 20 - 34cm, ensuring correct thigh support as the child grows.





Step 1



Step 2

// PATENTED QUICK-FOLD TECHNOLOGY

Folds into a very compact package in two simple steps (just pull!). The lightweight aluminium base weighs just 12kg and is easy to handle, fold and lock for safe transfers. This means less strain on your back.



Approved for use as an occupant's seat in a motor vehicle when used with an approved Wheelchair Tiedown and Occupant Restraint System in the front facing position.

EASY HANDLING FOR DAILY USE

Quickly attach and remove the seat from the stroller in one simple motion. This makes transport easy and ensures that the stroller fits in the smallest of cars.




// EASY TO CHANGE - FRONT AND REAR FACING

Face the seat rearwards to interact with your child or forward so they can explore their surroundings.

// SLIDE-N-LOCK SEAT ADAPTER

Removing and re-attaching the seat is simple with the Slide-n-Lock adapter. Rails guide the seat into position until an audible lock indicates that the seat has been safely attached.

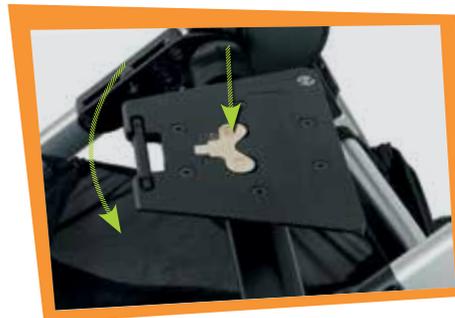



ZIPPIE voyage



COLOUR TO MATCH YOUR STYLE!

Choose one of four upholstery colours to complement the stylish mainstream Baby Jogger™ base.





FUNCTIONAL ACCESSORIES OF THE ZIPPIE VOYAGE

The Zippie Voyage is available with a large range of options and accessories designed to make daily tasks a breeze. For these accessories please contact your local Sunrise Medical Dealer.

Available through Sunrise Medical approved distributors



// UPPER EXTREMITY SUPPORT TRAY

This support tray is available in two versions – in light grey or with clear poly material. Both versions are depth and height adjustable.



// UTILITY UNDER FRAME TRAY

Versatile hard tray storage option for e.g. batteries or vents. Maximum load up to 23kg.



// FRONT MOUNT UTILITY PLATFORM

Front mount storage platform to enable easy access to carried objects.



// SHOE HOLDER

The shoe holders helps to improve overall positioning by holding the foot securely.

Available from Baby Jogger approved suppliers



// HOOK UTILITY

Useful hook to hold multiple bags. Easy to attach with Velcro strap that stays put and will not slide around.



// GLIDER BOARD

The glider board is a great way to let your older child hitch a ride. It attaches easily to the rear axle of your stroller and has a non-slip surface so you know your child will be safe.



// GRAB BAR

The removable grab bar gives your child a safe and comfortable area to hold on to whilst strolling. It can be angled forward or backwards.



// CUP HOLDER

For a water bottle or a coffee-to-go this cup holder uses a self-levelling mechanism to prevent unwanted spills when strolling.

TECHNICAL DATA

Zippie Voyage // Early intervention

Max. user weight:	34kg	Max. reversing width:	110cm
Seat width:	24 – 30cm, with insert down to 15.5cm (excl. upholstery)	Max. safe slope:	9° (with use of parking brake)
Seat depth:	20 – 34cm (excl. upholstery)	Rear wheel:	12.5" pneumatic with airless inserts
Lower leg length:	15.5 – 31cm	Front wheel:	7"
Backrest height:	41 – 61cm (excl. upholstery)	Seat weight:	6kg
Backrest angle range:	100° (80° to 180°)	Base weight:	12kg
Seat angle range:	35° (-5° to 30°)	Removable parts:	Seat and wheels - for transport only
Seat height:	59 – 66cm (3 steps possible 59 / 62 / 66cm)	Colour options:	4 colours - blue, red, beige or grey trims
Total width:	63cm	Frame material:	Aluminium
Total length:	110cm	Crash tested:	Yes; tested in accordance with Annex A of ISO 7176-19 (2008)
Total height:	98 – 105cm	Transport without user:	Secure the product safely in the vehicle (e.g. car or plane)
Folded length:	80cm	Intended use and environment:	For children with limited mobility: indoor and outdoor use, operated by caregiver
Folded height:	37cm		

For further information on the full specification, options and accessories, please refer to the order form.

All information is subject to change without notification. Please consult Sunrise Medical with any queries you may have.

If you are visually impaired, this document can be viewed in PDF format at www.SunriseMedical.co.uk or alternatively is available on request in large text.



Sunrise Medical

Thorns Road,
Brierley Hill,
West Midlands,
DY5 2LD

Tel.: +44 (0) 845 605 66 88

Fax: +44 (0) 845 605 66 89

E-Mail: enquiries@sunrisemedical.co.uk

www.SunriseMedical.co.uk